4-WHEEL SCOOTER




SUPERIOR COMFORT, PERFORMANCE AND STYLE

The **Celebrity® DX** delivers superb performance and ease of use with sleek, sporty styling for an overall package that impresses with elegant simplicity. Loaded with standard features like a deluxe seat with headrest and manual recline, and lower light for pathway illumination, the Celebrity DX is equipped to excite.





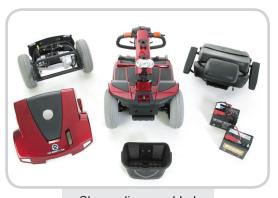

Shown disassembled

PRIDE® MOBILITY SCOOTERS
LIVE YOUR BEST®

SPECIFICATIONS

CELEBRITY® DX 4-wheel

Weight capacity:	159 kg (350 lbs.)
Maximum speed:	Up to 9.65 km/h (6 mph)
Ground clearance ² :	8.26 cm (3.25")
Turning radius:	149.86 cm (59")
Overall length ² :	118.75 cm (46.75")
Overall width ² :	61 cm (24")
Tyres:	Front: 9 cm x 25 cm (3.5" x 10") Rear: 10 cm x 25 cm (4" x 10") Anti-tips: Rear
Suspension:	Front & rear
Per-charge range ^{1,4} :	Up to 20.63 km (12.8 miles)
Total w/o batteries:	65.31 kg (144 lbs.)
Heaviest piece when disassembled:	Front section 29.02 kg (64 lbs.)
Standard Seating:	Type: Deluxe seat with headrest and manual recline Material: Black vinyl Dimensions: Width: 45.72 cm (18") Depth: 45.72 cm (18")
Drive System:	Rear-wheel drive, sealed transaxle, 24-volt DC motor
Brakes:	Regenerative and electromechanical
Battery requirements ³ :	(2) 12V deep cycle Size: 32 AH or 40 AH Weight ⁵ : 32 AH: 11.11 kg (24.5 lbs.) each 40 AH: 14.5 kg (32 lbs.) each
Battery charger:	Off-board, 5 amp

Colour:



1 Varies with user weight, terrain type, battery amp-hour (AH), battery charge, battery condition and tyre condition. These specification can be subject to a variance of (+/- 10%).

NOTE: All specifications subject to change without notice.

The information contained herein is correct at the time of publication; due to Pride's commitment to constant improvement and development, we reserve the right to alter specification without prior notice.

Theoretical range is calculated under test conditions in accordance with European standards. Actual range in normal use will depend on many factors, including the condition of the vehicle and its batteries, the weight of the driver, correct type pressures, ambient temperature, and the gradient and surface of the road or payment.

FEATURES

- Speed reduction technology automatically reduces travel speed while turning corners
- State-of-the-art rear suspension
- Top speed of up to 9.65 km/h (6 mph)
- Maximum weight capacity of 159 kg (350 lbs.)
- Ultra-simple, auto-latching lockup for quick disassembly and assembly
- Easily accessible tie down points (for transport of unoccupied scooter)
- Front basket and rearview mirror standard
- Tight turning radius of 149.86 cm (59")
- Includes standard angle-adjustable light
- Front frame-mounted seat post for optimal stability
- · High visibility reflector
- Microprocessor-based controller
- Tiller rotation lock permits tiller to be used as a handle during transport

OPTIONS

- · Cane/crutch holder
- Oxygen tank holder
- · Cup holder
- Rear basket
- Walker holder



Easy-drive tiller with wraparound handles







² Due to manufacturing tolerances and continual product improvement, this specificatio can be subject to a variance of (+ or - 3%).

³ AGM or gel cell type required.

⁴ Tested in accordance with ANSI/RESNA, WC Vol2, section 4 & ISO 7176-4 standards. Results derived from theoretical calculations based on battery specification and drive system performance test conducted at maximum weight capacity

⁵ Battery weight may vary based on manufacturer.