




SKU: S11291

AMBITURN SIT TO STAND TRANSFER AID

\$1,299.99

The Ambiturn Sit to Stand Transfer Aid assists with pivoting transfers between different seat surfaces and can transport users' short distances. The product requires users to be weight bearing and able to push themselves into a standing position. Includes a low-profile base which provides substantial stability with rear bar for attendees to counterbalance. The Ambiturn has handles with multiple gripping points for varying support and comfort, and padded lower leg supports which can be adjusted.

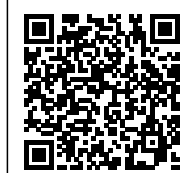

Scan QR code to
view product online

PRODUCT DESCRIPTION

Key Features

- Requires user to pull on handle to stand up and to be weight bearing
- Narrow, recessed front fits between wheelchair wheels and around toilet for perfect positioning
- Low profile base facilitates user feet placement
- Substantial base area provides stability, with rear bar also allowing carer to counterbalance
- Handle has multiple gripping points for varying carer/user support requirements and is angled to enable “nose over toe” standing movement
- Opening on central handle bar allows quick fitting of optional support belts
- 4 front and rear castors combined with 2 x 13cm (5”) central wheels ensure effortless manoeuvring even with heavier users
- Foot operated brakes on centre wheels provide stability during standing/sitting movement
- Padded lower leg supports can be adjusted (tool free) in width, angle and height to suit user comfort
- Easily dismantled/reassembled without tools for transportation and storage
- Easy wipe clean surfaces

ADDITIONAL INFORMATION

Brand	Drive Medical
Overall Width	57.5cm
Overall Depth	63cm
Overall Height	118cm
Base Plate Length	63cm
Base Plate Width	57.5cm
Leg Support Height Range	30cm – 47cm
Unit Weight	16.5Kg
Weight Capacity/Max SWL	150kg