



**SKU:** S11499

## **BELLAVITA BATH LIFT**


**\$999.99**

**EXCLUSIVE TO ILS.**

**Enjoy a bath safely and comfortably with the #1 rated bath lift in UK, the Drive Medical Bellavita reclining bath lift.**

Made in Germany, the Bellavita bath lift is the lightest battery-powered reclining bath lift available, weighing only 9.3kg. Developed with the help of nursing professionals, it is ergonomically designed for you to enjoy a great soak in the bath with a record seat height of 48cm, very low seat travel of 6cm and a backrest that reclines to a market leading 50 degrees.

Built with many safety features, this ultra-modern bath lift is much loved by our customers as a daily living aid.

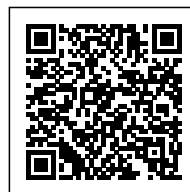

Scan QR code to  
view product online

## PRODUCT DESCRIPTION

Enjoy a bath with confidence using the lightest battery powered bath lift available - the Bellavita bath lift. At the touch of a button using the hand control of this ergonomically-designed bath lift, users with restricted mobility are gently lowered down into the bath. The bath lift also lifts the user to the top again for safely getting out.

### **Key features include:**

- At the push of a button, it lowers persons with limited mobility to the bottom of the tub for bathing and relaxation, then elevates them back to the top of the tub for a safe exit.
- The backrest reclines to a market-leading 50 degrees to accommodate customers requiring a greater hip-to-back sitting position, such as those with stiff hips.
- The cutout in the seat's front facilitates personal hygiene.
- Can be used with or without covers; washable hygienic covers are simply fitted/removed and offer additional bathing comfort (available in white or blue)
- Extremely low seat travel of 6 cm (2.4") from the bottom of the bath enhances water immersion and is perfect for shallow baths.
- Record seat height of 48 cm (18.9") allows for safe level transfers in deeper baths
- Weighs only 9.3 kg (20.5lb) to make transporting and handling easier.
- Quickly and simply separates into two lightweight halves for simple installation/removal, storage, and travel
- Contoured, slim-line profile facilitates handling and provides a contemporary, aesthetically pleasant appearance that complements most bathrooms.
- Four quick-release suction feet permit combat-free removal from bath
- Lightweight, completely waterproof, floatable hand control with easy-to-use buttons
- The handheld remote's lithium-battery technology allows for convenient recharging in numerous locations.
- An indicator light illuminates when charging is necessary, and a safety mechanism stops the bath lift from lowering without sufficient energy to return it to the top of the bath.
- The emergency stop button built into the handle ensures maximum safety
- User-friendly jack plug connection to charger for simplicity of usage
- Produced from high-quality, lightweight plastic that is recyclable and simple to clean
- Extensive array of available accessories to optimise bath installation and serve a variety of consumers
- Weight limit of 140 kg

**Additional item can be purchased separately:**

- [Washable hygienic covers](#) – Available in blue, the covers are easily fitted/removed and provide extra comfort when bathing.

## ADDITIONAL INFORMATION

<b>Brand</b>	Drive Medical
<b>Max user weight (SWL)</b>	140kg
<b>Total weight (excl. hand control)</b>	9.2kg
<b>Overall length</b>	83cm (32.5 inches) – with backrest in a reclined position
<b>Overall width</b>	70cm (27.5 inches) – seat including side flaps
<b>Seat width</b>	35cm (14 inches)
<b>Seat depth</b>	42cm (16.5 inches)
<b>Seat length</b>	57cm (22.5 inches)
<b>Seat height</b>	6cm – 48cm (2.5 – 19 inches)
<b>Backrest height</b>	64cm (25 inches)
<b>Backrest width</b>	32 – 36cm (13 – 14 inches)
<b>Backrest reclines</b>	50 degrees
<b>Baseplate width</b>	31cm (12 inches)
<b>Baseplate length</b>	56cm (22 inches)
<b>Warranty</b>	2 years warranty on battery/handset.
<b>Colour</b>	Blue