


SKU:

## KCARE COMPACT ELECTRIC RECLINER LIFT CHAIR

~~\$4,199.99~~ \$4,199.99

Designed for customers who require assistance sitting and leaving their chair, the K Care Compact Lift Chair offers proven pressure relief utilising the unique Air Comfort Seating System. Removable seat and back cushions ensure a snug fit while the electric lift and recline mechanism aids easy operation. Features VersoTela Fusion fabric is designed specifically for pressure care and sensitive skin

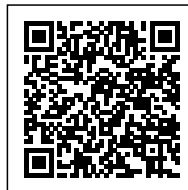

Scan QR code to  
view product online



## VARIATIONS

Image	SKU	Price	Description	Sizes
	S2430	\$3,999.99	Single Motor Small	
	S1469	\$3,999.99	Single Motor Medium	
	S7527	\$3,999.99	Single Motor Large	
	S2236	\$4,199.99	Twin Motor Small	

Image	SKU	Price	Description	Sizes
	S1039	\$4,199.99	Twin Motor Medium	
	S3190	\$4,199.99	Twin Motor Large	

## PRODUCT DESCRIPTION

The Air Comfort seating system is designed to aid in the treatment and prevention of pressure problems related to long term seating. Utilising 4 pressure-less, air filled, high frequency welded PVC bags, Air Comfort products work to distribute the occupants weight, thereby ensuring there are no main pressure points restricting circulation.

The Air Comfort Kcare lift chair comes with front and lockable rear castors, allowing easy maneuvering when required.

### Features

- Total surface air pressure system
- Removable seat and back cushions
- Electric lift and recline mechanism
- Front and rear castors
- Rear locking castor
- Wall Saver mechanism
- Wooden hand knuckles
- Back-up battery capability
- Proven pressure relief for long-term seating
- Assisted sitting, lifting, and reclining
- Easy side patient transfer
- Easy to operate and manoeuvre
- Easy to clean
- Easy to repair and maintain

### Single vs twin vs quad motor lift chairs - what's the difference?

All ILS lift chairs have a lifting action, to support the user to stand from a seated position.

- Single motor chairs move the back rest and leg rest at the same time.
- Twin motor chairs allow the backrest and leg rest to move independently of each other, typically allowing for a greater degree of leg elevation and greater degree of recline compared to a standard single motor chair.
- Quad motor chairs allow independent movement of headrest, lumbar, backrest and footrest allowing for the greatest level of re-positioning of all chairs.



## ADDITIONAL INFORMATION

<b>Brand</b>	KCare
<b>Max SWL</b>	140 kg
<b>Back Height</b>	Small: 680 mm / Medium: 740 mm / Large: 800 mm
<b>Seat Height</b>	Small: 480 mm / Medium: 520 mm / Large: 520 mm
<b>Seat Depth</b>	Small: 440 mm / Medium: 510 mm / Large: 510 mm
<b>Seat Width</b>	Small: 500 mm / Medium: 500 mm / Large: 600 mm
<b>Overall Width</b>	Small: 750 mm / Medium: 750 mm / Large: 820mm
<b>Overall Height</b>	1060 mm
<b>Manufacturer's part number</b>	AC59018
<b>Sizes</b>	Single Motor Small, Single Motor Medium, Single Motor Large, Twin Motor Small, Twin Motor Medium, Twin Motor Large
<b>Warranty</b>	A warranty of 5 years on the frame, 2 years on castors and motor 1 year on fabric and struts.