




SKU:

## **FORTE SOVEREIGN PLUS MATTRESS (ILS)**

**~~\$1,499.99~~ \$1,499.99**

The Forte Sovereign Plus Mattress is a medium-risk static support surface. Made in Australia, it combines visco-elastic immersion with zoned comfort and a durable cover for use across home, hospital and aged care environments.

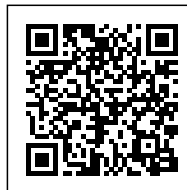

Scan QR code to  
view product online

VARIATIONS

Image	SKU	Price	Description	Select Size
	S16394	\$1,349.99	Single	
	S16395	\$1,349.99	Long Single	
	S16397	\$1,499.99	King Single	
	S16396	\$1,499.99	Long King Single	

## PRODUCT DESCRIPTION

The Sovereign Plus Mattress is engineered to provide effective support and comfort for those at medium risk of pressure injury. Its visco-elastic upper layer conforms to pressure points, while the dual-density secondary layer enhances immersion and envelopment.

Profiled stabilising sides assist with positioning and bed conformity. A firm base provides sacral support and assists with repositioning.

The moisture-resistant Premiflex cover is stretchable, anti-shearing and easy to clean – making it suitable for high-demand settings.

### Who it's for

- For clients with medium pressure injury risk who need enhanced support and hygiene.
- For therapists requiring a durable, clinically supportive mattress with versatile use across facilities or at home.

### Key features

- **Visco-elastic support surface** - sensitive immersion interface, dual density secondary immersion support layer, zoned for optimal immersion and comfort
- **Profiled stabilising sides** - enables ease of primary patient care, facilitates patient transfers, unique hinge points allow the mattress to easily conform to varied bed positions, encourages central patient positioning
- **Firm base support layer** - provides sacral protection and exceptional comfort and support for repositioning and patient care
- **Sloped heel relief** - enhanced pressure redistribution and protection for critical at-risk heel region
- **Premiflex cover** - superior infection control, hygiene and high-performance cover that is suitable for demanding care settings. Produced in Belgium, the cover is a 4-way stretch PU interface which is moisture vapour permeable and anti-shearing. Heavy-duty upgraded PVC base material protects premature cover deterioration

