




SKU: S18675

INVACARE BIRDIE EVO 180KG PATIENT LIFT HOISTS

\$4,699.99

The Invacare Birdie Evo 180 kg Patient Lift Hoist combines ergonomic design with advanced lifting features to support comfortable and dignified transfers in all care settings. It offers tool-free folding, 360° rotation and electric leg adjustment for flexibility and ease of use.

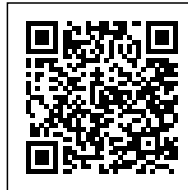

Scan QR code to
view product online

PRODUCT DESCRIPTION

The Invacare Birdie Evo 180 kg Patient Lift Hoist is engineered to deliver maximum safety and comfort during everyday lifting tasks. Whether transferring from a bed, chair or the floor, the Birdie Evo offers smooth, quiet operation with thoughtful features that support both the person's dignity and carer efficiency.

The high boom design provides generous space around the person, while the front clearance reduces the risk of contact with the knees. Its electric leg adjustment allows for positioning flexibility between furniture, with a spread range of 64-102 cm.

The full 360° rotation allows for smoother, more controlled positioning. An improved ergonomic hook design makes sling attachment faster and more secure, while the Smartlock system allows quick, tool-free spreader bar changes for added convenience.

The integrated Slow'R dampener reduces swaying during lifts, improving comfort and control. The unit folds without tools and features a detachable battery for simple maintenance.

Who it's for

- For clients who require full-body hoisting with enhanced positioning and comfort.
- For therapists wanting a versatile, ergonomic hoist for safe use across hospital, care and home settings.

Key features

- **360° rotation** – simplifies positioning and reduces anxiety during transfers
- **Electric leg adjustment** – accommodates varying furniture widths
- **High boom clearance** – increases comfort and reduces contact with knees
- **Tool-free folding** – enables quick and compact storage
- **Improved sling hook design** – secure, one-handed attachment with reduced risk of sling detachment
- **Smartlock spreader bar system** – quick-release mechanism for fast, safe bar changes
- **Integrated Slow'R dampener** – minimises rocking motion while in use
- **Ergonomic push bar** – enhances manoeuvrability for the carer
- **Detachable battery pack** – simplifies maintenance
- **Includes 2-point spreader bar** – compatible with optional 4-point accessory (S300789)

ADDITIONAL INFORMATION

Brand	Invacare
Minimum Base Width	640 mm
Maximum Inside Base Width (Max Spread)	1020 mm
Lifting Range	415 mm-1660 mm
Length	1220 mm
Underbed Clearance	110 mm
Product Weight	40 kg
Battery Time	40 Cycles
Maximum Weight Capacity	180 kg