




SKU: S18674

INVACARE BIRDIE COMPACT PATIENT LIFT HOIST

From: \$3,699.99

The Birdie Compact is only 55cm wide making it the only lifter that is TÜV certified and is this narrow, a unique lifter to the market. The ergonomically designed push bar allows the carer to position their hands comfortably when moving the lifter, reducing the likelihood of strain. The 100mm smoothing rolling castors make it easy to push and pull the lifter, reducing the force required to manoeuvre.

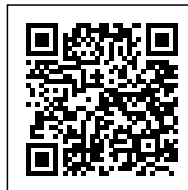

Scan QR code to
view product online

PRODUCT DESCRIPTION

- Birdie Evo Compact is simple to fold, unfold, and dismantle without tools. This makes it very easy to handle. The hoist takes up minimal space when stored and is easy to push and transport. If required, the mast and boom can also be easily detached from the base for easy transportation without any tools required. Below are some key features of the Invacare Birdie Compact:
- Emergency stop, lowering, and integrated charger: Manual or electric emergency lowering (compact version only electric lowering). The control box is supplied with a detachable or fixed battery.
- Manual leg spread: Users can simply step on the pedal to adjust the leg spread.
- Long boom – high user comfort: The client can easily be rotated and positioned. This reduces the risk of injury to the knees.
- Easy to fold and dismantle: Allows for ease of transportation and storage. No tools required when folding/unfolding or when dismantling the hoist.
- Mounted spreader bar: This Allows correct and safe positioning of the spreader bar.

ADDITIONAL INFORMATION

Brand	Invacare
Minimum Base Width	520 mm
Maximum Inside Base Width (Max Spread)	815 mm
Lifting Range	530 mm-1590 mm
Length	1080 mm
Underbed Clearance	110 mm
Product Weight	31 kg
Battery Time	40 Cycles
Maximum Weight Capacity	150 kg