




SKU:

ROYALE MAYFAIR SELECT VINYL (SINGLE MOTOR)

\$1,399.99

One of our most popular lift chairs is the affordable Royale Medical's Mayfair Select electric recliner lift chair. It is much loved by our customers for its quiet and smooth three position single motor lift system and very comfortable 3-section *waterfall* backrest design for optimum comfort and back support. Available in four classic colours, it is made with select studio high-quality vinyl to suit any setting or home decor. Enjoy the comfier way of life with the Mayfair Select electric recliner lift chair!

Note: Armrest covers are sold separately.

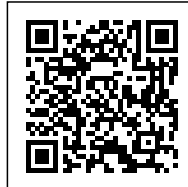

Scan QR code to
view product online

VARIATIONS

Image	SKU	Price	Description	Colours
	S12203	\$1,399.99		Chocolate Brown
	S12204	\$1,399.99		Charcoal Black
	S12205	\$1,399.99		Vanilla Beige
	S12206	\$1,399.99		Mushroom Brown

PRODUCT DESCRIPTION

Royale Medical's Mayfair Select electric recliner lift chair offers support, comfort and style for every day use at an affordable price.

Key features include:

- Very comfortable 3-section *waterfall* backrest design for optimum comfort and upper, middle and lower back support.
- Easy-to-use functions so you can easily recline the backrest and legrest.
- Quiet and smooth single motor lift system.
- Quality actuators and American-made mechanisms allows motors to run quietly for a smooth and discreet operation.
- Quality foams and no sag springs, providing you with supreme postural support.
- High-quality vinyl upholstery in four colours – Choose from charcoal black, vanilla beige, mushroom brown, chocolate brown.
- Backup power supply for a peace of mind.
- Manufacturer's warranty period of 2 years for electrical and steel lift/recliner mechanism.

Single vs twin vs quad motor lift chairs - what's the difference?

All ILS lift chairs have a lifting action, to support the user to stand from a seated position.

- Single motor - chairs move the back rest and leg rest at the same time.
- Twin motor - chairs allow the backrest and leg rest to move independently of each other, typically allowing for a greater degree of leg elevation and greater degree of recline compared to a standard single motor chair.
- Quad motor - chairs allow independent movement of headrest, lumbar, backrest and footrest allowing for the greatest level of re-positioning of all chairs.

ADDITIONAL INFORMATION

Brand	Royale Medical
Max SWL	150 kg (330.6 lbs)
Motor	Single
Overall height	105 cm (41.3 inches)
Overall width	75 cm (29.5 inches)
Overall depth	75 cm (29.5 inches)
Armrest height (from floor)	62 cm (24.4 inches)
Seat height (from floor)	48 cm (18.8 inches)
Seat width (measured on sitting area)	52 cm (20.4 inches)
Seat depth (measured on sitting area)	50 cm (19.6 inches)
Seat max height (lifted)	56 cm (22 inches)
Backrest length (fully reclined)	32 cm (12.5 inches)
Max legrest height (lifted)	46 cm (18.1 inches)
Foam density	28/70
Upholstery	Studio vinyl
Colour options	Charcoal black, vanilla beige, mushroom brown, chocolate brown
Warranty	2 years manufacturer's warranty on the frame integrity, workmanship and electronics.