

SKU: S18612



MINILIFT200 PATIENT LIFTER HOIST

\$6,995.99

MiniLift200 is a mobile sit-to-stand lift which has been developed to, as gently as possible, assist the user when rising from a sitting to a standing position. When MiniLift200 is combined with the appropriate lifting accessories, the user gets support under the feet, for the front of the lower legs and behind the back, which provides for a safe and secure, yet active sit-to-stand procedure.

To be able to understand how comfortable MiniLift200 is compared with other sit-to-stand devices on the market, it must be tried! The construction safely moves the user forwards and upwards in a natural movement pattern and, at the same time, the leg muscles and balance are exercised.


MiniLift200 is manufactured in steel for stability and strength and is suitable also for tall and heavy users weighing up to 200 kg/440 lbs.

Features:

- Unique solution – extremely comfortable, active sit-to-stand with natural movement patterns

- Small and convenient
- Easy to manoeuvre
- Ergonomic and easy-to-use hand control with 4 buttons
- Quick guide on the lift
- Grip-friendly handles with multiple options for hand placement
- Soft, adjustable lower leg support
- Low foot plate with anti-slip surface
- Easy rolling, maintenance-free castors
- Lockable rear castors
- High capacity battery with long service life
- Control box with diagnostics, e.g. lift counter and overload indicator
- Electrical emergency stop readily available on the control box

MiniLift200 meets the requirements of the Medical Devices Directive for Class 1 products 93/42/EEC.

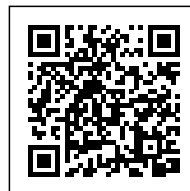

Scan QR code to
view product online

ADDITIONAL INFORMATION

Brand	SystemRoMedic
Safe working load (SWL)	200 kg / 440 lbs
User length, min-max	1400-2000 mm / 55.1-78.7"
Material	Powder-coated steel
Height	1000 mm / 39.4"
Length	980 mm / 38.6"
Width of base, outer measurement	590-1030 mm / 23.2-40.6"
Width of base, inner measurement	450-875 mm / 17.7-34.4"
Weight	49 kg / 108 lbs
Base widening	Electrical
Height of footplate, for user	95 mm / 3.7"
Castor size, diameter	100 mm / 3.9"
Lifting speed, without load	35 mm / 1.4" per sec
Emergency lowering	Manual and electrical