




SKU:

## JAY ION CUSHION

~~\$649.99~~ **\$649.99**

The JAY Ion is a comfortable skin protection and positioning wheelchair cushion composed of a multi-layered foam base, moisture resistant inner cover, and an anti-microbial outer cover. The top layer of soft, Visco memory foam featured in this wheelchair cushion conforms to each individual and distributes weight over a larger area for pressure reduction and maximum comfort.

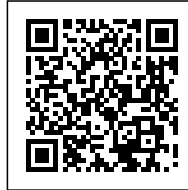

Scan QR code to  
view product online

## VARIATIONS

Image	SKU	Price	Description	Size
	S16031	\$549.99		14x15"
	S17397	\$549.99		15x15"
	S16525	\$549.99		15x16"
	S17145	\$549.99		15x18"

Image	SKU	Price	Description	Size
	S12670	\$549.99		16x16"
	S14928	\$549.99		16x17"
	S15092	\$549.99		16x18"
	S17301	\$549.99		16x19"

Image	SKU	Price	Description	Size
	S13270	\$549.99		17x17"
	S13900	\$549.99		17x19"
	S2410	\$599.99		18x16"
	S7729	\$599.99		18x17"

Image	SKU	Price	Description	Size
	S12237	\$599.99		18x18"
	S11264	\$599.99		18x19"
	S16999	\$599.99		18x20"
	S14397	\$599.99		18x22"

Image	SKU	Price	Description	Size
	S14531	\$649.99		20x18"
	S12038	\$649.99		20x20"

## **PRODUCT DESCRIPTION**

Dual layered foam base with varying firmness

- The Ion wheelchair seat cushion features a dual-layered foam base for a combination of structural support and comfort. For optimised immersion for each individual, the foam's firmness is determined by the cushion width, which directly corresponds with client weight.

Visco memory foam top layer

- The top layer of soft, Visco memory foam featured in this wheelchair cushion conforms to each individual and distributes weight over a larger area for pressure reduction and maximum comfort.

Innovative dual cover system

- For convenience and protection, the JAY Ion features a dual-cover system. The inner cover is water-resistant and easy to wipe off and clean, while its Aqua-guard zipper and anti-wicking thread protect the foam base from fluids. The machine-washable, stretch outer cover reduces surface tension and utilises naturally anti-microbial, silver-impregnated Ionic+ A layer of 3DX spacer fabric promotes airflow and dissipates heat and moisture for additional skin protection.

Common client characteristics

- May be at moderate risk of skin breakdown
- May be able to perform independent weight shifts
- May have moderate positioning needs
- May be in a wheelchair for extended periods of time

## **ADDITIONAL INFORMATION**

<b>Max SWL</b>	136 kg : 14" -21" wide, 225 kg: 22" – 24" wide
<b>Product Weight Max</b>	1.5 kg(based on 18" x 18" size)
<b>Cover Type</b>	X-Static® with Ballistic Nylon
<b>Brand</b>	Sunrise Medical