



SKU:


## ROYALE MULTI-MOTION BED - KING SINGLE

**\$3,699.99**

Introducing the Royale Multi-Motion adjustable home care bed, carefully crafted to meet the clinical needs of users while seamlessly integrating into any home environment. With full hospital-grade profiling functions, it prioritises user comfort, positioning, and safety.

Available in various sizes with a range of accessories to cater to an individual's functional and clinical requirements. In addition to its functionality, this bed features a modern design that effortlessly complements any home decor, ensuring seamless integration into your living space.

Available in 4 sizes to suit the size of your living space: Single, King Single, Double, and Queen.

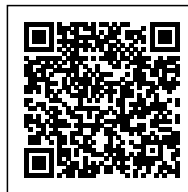

Scan QR code to  
view product online

## VARIATIONS

Image	SKU	Price	Description	Colour
	S306613	\$3,699.99		Black
	S306614	\$3,699.99		Grey

## PRODUCT DESCRIPTION

The movement features of the Royale Multi-Motion Bed, are designed to enhance comfort, safety, and convenience.

- **Height:** Raise and lower the entire bed to enhance accessibility for the user and carer, providing better ergonomics during user care while reducing the risk of poor manual handling.
- **Knee Break:** Allows independent elevation of the bed's knee area. Helps improve circulation and alleviate lower back and leg pressure, providing personalised comfort for various activities like reading, watching TV, or relaxing.
- **Back Rest:** Enables independent adjustment of the upper bed section, including the head and upper body, for personalised positioning. Enhances comfort, postural support, and ease of transfers.
- **Trendelenburg:** Allows a user to be positioned in supine (lying flat on their back) with their feet elevated above their head, or in reverse Trendelenburg, with their head elevated higher than their feet.

Additional features:

**Illuminated Hand control:** Hand control features illuminated buttons and intuitive icons to allow for easy adjustments, even in low-light conditions.

**Heavy duty castors:** Adjust positioning with ease. Equipped with eight smooth-swivel castors, with one lockable castor in each corner.

**Upholstery:** Industrial grade, soft-touch woven fabric is available in 2 contemporary colours: Ardo Grey and Jet Black to suit your personal style.

Available Accessories:

A range of easy to attach accessories are available to provide additional safety and support:

- [Headboard](#)
- [Footboard](#)
- [Wide bedside rail](#)
- [Three quarter length bedside rail](#)
- [Bed stick](#)
- [IV Pole](#)

