




**SKU:** S2229

## **SHOWER COMMODE - MOBILE - SELF PROPELLED - FOLDING**

**From: \$2,499.99**

The folding shower commode assists with the transportation of a patient for showering and toileting. It is easily positioned over the toilet for toileting and folds neatly for transportation or storage where space is limited. The padded commode seat (available as open or closed front) features adjustable height swing-away leg-rests assist patient transfers as do it's swing-back arms for side transfer. The hinge point is in the middle of the seat base of the shower commode allowing it to fold inwards. Australian made.

**Seat sold separately**

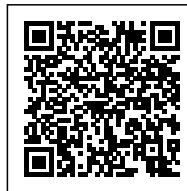

Scan QR code to  
view product online

## **PRODUCT DESCRIPTION**

A removable PVC slatted backrest aids in storage and patient access and the folding shower commode moves on a pair 2 x 125 mm castors and 2 x 22" rear wheels. Built with a durable stainless steel frame. Padded Seat is optional.

## **ADDITIONAL INFORMATION**

<b>Seat Width</b>	445 mm
<b>Seat Depth</b>	430 mm
<b>Seat Height</b>	590 mm
<b>Arm Height</b>	240 mm
<b>Back Height</b>	365 mm
<b>Seat to Footplate</b>	380 mm – 520 mm
<b>Overall Width</b>	610 mm
<b>Max SWL</b>	120 kg
<b>daily-living</b>	Commodes