




SKU:

## QUATTRO ACUTE MATTRESS (FORMERLY TALLEY)



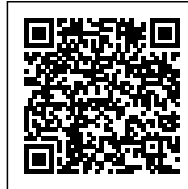

**~~\$14,995.99~~ \$14,995.99**

The Quattro Acute mattress system is a market leading and highly regarded 1-in-4 alternating therapy mattress. The QUATTRO Acute is unique in offering DEEP CELL THERAPY and this top of the range support surface is ideal for your most dependent patients. Designed with TISSUEgard, Ortho-Differential Support and a multi-stretch cover to provide patients with pressure relief to prevent and manage patients at risk of or with existing pressure injuries.



Scan QR code to  
view product online

## VARIATIONS

| Image  | SKU     | Price       | Description | Size and Cells            |
|--|---------|-------------|-------------|---------------------------|
|  | S306784 | \$12,995.99 |             | Single 27 cells           |
|  | S311257 | \$14,995.99 |             | King Single 27 cells      |
|  | S306783 | \$12,995.99 |             | Single Long 28 cells      |
|  | S311259 | \$14,995.99 |             | King Single Long 28 cells |

## **PRODUCT DESCRIPTION**

Quattro Therapy Systems provide a clinical preventative measure to provide pressure relief and manage patients at risk of developing or with existing pressure Injuries. Please see below some of the product's key features

- DEEP CELL THERAPY allows the cells of the QUATTRO Acute to run at lower internal pressures. This minimises the pressure applied to the patients skin and subcutaneous tissues during the cycle.
- The active 1-in-4 cell cycle supports 75% of the patient's body at all times, on groups of three inflated cells, whilst the fourth cell deflates sufficiently to redistribute the pressure and encourage tissue reperfusion. Additional benefits include enhanced patient comfort, less awareness of support surface movement and a reduction of stimulus-related complications such as muscle spasm.
- TISSUEgard enables the partial immersion and envelopment of the patient into the support surface, reducing the pressure on the patient's skin and decreasing the pressure differential between inflated and deflated cells, which helps reduce tissue strain and associated shearing forces.
- ODS results in firmer outer edges of the mattress which facilitate patient transfers and provide extra support, safety and comfort for larger or heavier patients. It also creates a softer central area of the mattress, ideal for smaller, lighter patients.

## **Features**

- Active 1-in-4 cell cycle
- Deep Cell Therapy
- TissueGard
- Ortho-differential Support
- Multi-stretch cover

- Microclimate management
- Fully launderable
- Simply designed interface
- Fast, safe and easy manual interventions
- Continuous low-pressure mode
- Max inflation
- Adjustable comfort control

## ADDITIONAL INFORMATION

|                                     |                         |
|-------------------------------------|-------------------------|
| <b>Brand</b>                        | Direct Healthcare Group |
| <b>Power Unit Width</b>             | 346 mm                  |
| <b>Power Unit Height</b>            | 258 mm                  |
| <b>Power Unit Depth</b>             | 156 mm                  |
| <b>Power Unit Weight</b>            | 4 kg                    |
| <b>Power Unit Cycle Time</b>        | 16 mins                 |
| <b>Mattress Length</b>              | 1980 mm                 |
| <b>Mattress Width</b>               | 900 mm                  |
| <b>Mattress Depth</b>               | 240 mm                  |
| <b>Mattress Weight</b>              | 12.5 kg                 |
| <b>Inner Construction Materials</b> | PU air cells            |
| <b>Cover Material</b>               | PU coated stretch nylon |
| <b>Number of Cells</b>              | 27, 28                  |
| <b>Max. User Weight</b>             | 250 kg                  |
| <b>Min. User Weight</b>             | None                    |