




SKU:

QUATTRO PLUS MATTRESS (FORMERLY TALLEY)

~~\$9,999.99~~ \$9,999.99

The Quattro Plus mattress system is a market leading and highly regarded 1-in-4 alternating therapy mattress. Designed with TISSUEgard, Ortho-Differential Support and a multi-stretch cover to provide patients with pressure relief to prevent and manage patients at risk of or with existing pressure injuries.

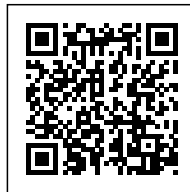

Scan QR code to
view product online

VARIATIONS

| Image | SKU | Price | Description | Size and Cells |
|--|---------|------------|------------------------------|---------------------------|
|  | S311260 | \$9,999.99 | | King Single Long 32 cells |
|  | S311256 | \$9,999.99 | | King Single 30 cells |
|  | S306785 | \$7,999.99 | | Single Long 32 cells |
|  | S306786 | \$7,999.99 | Contact us for availability. | Single 30 cells |

PRODUCT DESCRIPTION

Quattro Therapy Systems provide a clinical preventative measure to provide pressure relief and manage patients at risk of developing or with existing pressure Injuries. Please see below some of the product's key features.

- The active 1-in-4 cell cycle supports 75% of the patient's body at all times, on groups of three inflated cells, whilst the fourth cell deflates sufficiently to redistribute the pressure and encourage tissue reperfusion. Additional benefits of the 1-in-4 system may include increased patient comfort, and less awareness of the support surface movement.
- TISSUEgard™ enables the partial immersion and envelopment of the patient into the support surface, reducing the pressure on the patient's skin and decreasing the pressure differential between inflated and deflated cells, helping to reduce associated shearing forces.
- Ortho-Differential Support provides firmer outer edges of the mattress to facilitate patient transfers and provide extra support, whilst creating a softer central area of the mattress
- The multi-stretch cover promotes tissue offloading during cell deflation and reduces shear and friction during patient movement or manual repositioning. The material is waterproof, moisture-vapour permeable and all seams are welded to protect the inside of the mattress from fluid ingress.

Features

- Active 1-in-4 cell cycle
- Multi-stretch cover
- TISSUEGard
- Ortho-differential support
- Microclimate management
- Fully launderable
- Simply designed interface
- Fast, safe and easy manual interventions

- Patient Transport Facility
- Continuous low pressure
- Max inflation
- Adjustable comfort control

ADDITIONAL INFORMATION

| | |
|-------------------------------------|--|
| Brand | Direct Healthcare Group |
| Power Unit Width | 346 mm |
| Power Unit Height | 258 mm |
| Power Unit Depth | 156 mm |
| Power Unit Weight | 4 kg |
| Power Unit Cycle Time | 16 mins |
| Mattress Length | 1950 mm |
| Mattress Width | 880 mm |
| Mattress Depth | 180 mm |
| Weight | 9.4 kg |
| Inner Construction Materials | PU air cells |
| Cover Material | PU coated stretch nylon |
| Number of Cells | 30, 32 |
| Max. User Weight | 200 kg |
| Min. User Weight | None |
| Size | Single, Single Long, King Single, King Single Long |