



SKU: S2298


RAILS - TOILET RAILS THRONE 3- IN-1 STANDARD

\$599.99

Designed for safety, comfort, and benefits from the average build to the larger build individual.

This Rail enables the user to have support close to the body, where it is most needed. This model has a weight loading of 150 kg. The user can lower safely onto the toilet, and pushing up from the seat is also much easier. The Rail is safer than over frames as it is rigidly fixed to the toilet pan and does not need to be removed for other users.

Handles can also be folded down to facilitate a wheelchair transfer. Additionally, either one or both handles can be removed to suit some requirements.

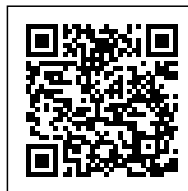

Scan QR code to
view product online

PRODUCT DESCRIPTION

The Standard 3-in-1 throne rail system is designed for a wide range of users. Providing excellent additional support for users when raising and lowering. Compatible with the entire range of throne accessory spacers. Arms can be positioned in a standard position, splayed position or folded down to provide optimal flexibility for varying individual needs and preferences.

ADDITIONAL INFORMATION

Width Between Handles	510 mm or 640 mm
Handle Height	260 mm
Standard Back Plate Internal	370 mm