




**SKU:** S7914

## **DRIVE MEDICAL TITAN POWER CHAIR**

**\$2,999.99**

The Titan is also great indoors, navigating tight, enclosed spaces with ease.

- Designed for indoor and outdoor use  
Front wheel drive system for optimum manoeuvrability
- Detachable padded and adjustable armrests to suit the comfort of the user
- Adjustable length controller mount
- Folding backrest for easy storage and transportation
- Height adjustable and reclining seat to suit comfort preferences
- Swivel seat facilitates easy transfers on and off the chair
- Height adjustable and removable headrest to suit user requirements
- The de-clutch operation to enable you to safely freewheel the chair where and when required
- Padded upholstery for additional comfort and support

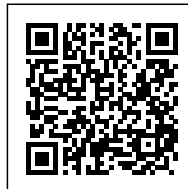

Scan QR code to  
view product online

## **PRODUCT DESCRIPTION**

The front wheel drive Titan Power chair is a comfortable and stylish option for users seeking a transportable power mobility product that's ideal for indoor and outdoor use. The Titan's front-wheel system positions the drive wheel forward of the user's centre of gravity, making the chair extremely stable and giving it great performance over small obstacles.

## ADDITIONAL INFORMATION

<b>Brand</b>	Drive Medical
<b>Width</b>	63.5cm
<b>Length</b>	102cm
<b>Height</b>	112cm
<b>Seat Width</b>	46cm
<b>Seat Depth</b>	46cm
<b>Seat Height</b>	55cm
<b>Weight</b>	72kg
<b>Max SWL (Maximum Weight)</b>	136kg
<b>Heaviest Piece Disassembled</b>	25 kgs
<b>Seat Weight</b>	19.35 kgs
<b>Battery Weight (PAIR)</b>	19.9 kgs.
<b>Base Weight</b>	25kg's
<b>Charger</b>	3A Offboard
<b>Maximum speed</b>	6.48 km/h
<b>Maximum range</b>	24 km
<b>Width between the armrests</b>	43cm - 53cm
<b>Warranty</b>	2 years