

SKU: S2438




TOPFORM ASHLEY ELECTRIC RECLINER LIFT CHAIR MEDIUM

\$2,499.99

The Topform Ashley Lift Chair is a luxury recliner chair designed with great features including a Tasmanian Oak frame with four optional (37mm) castors, high density foam padded seat and backrest with 3 separate padded horizontal cushions. Australian made. Wide variety of colour/fabric choices available.

For colour choices please email
onlineorders@ilsau.com.au

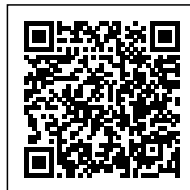

Scan QR code to
view product online

VARIATIONS

Image	SKU	Price	Description	Motor
	S2438	\$2,499.99		Single

PRODUCT DESCRIPTION

Full length, fully padded and upholstered armrests and available in a standard range of fabrics, non-standard fabrics at an additional cost. A two year frame, motor, transformer and hand control warranty and 12 month springs and component warranty are standard across the Topform Lift Chair range.

Single vs twin vs quad motor lift chairs - what's the difference?

All ILS lift chairs have a lifting action, to support the user to stand from a seated position.

- Single motor chairs move the back rest and leg rest at the same time.
- Twin motor chairs allow the backrest and leg rest to move independently of each other, typically allowing for a greater degree of leg elevation and greater degree of recline compared to a standard single motor chair.
- Quad motor chairs allow independent movement of headrest, lumbar, backrest and footrest allowing for the greatest level of re-positioning of all chairs.

ADDITIONAL INFORMATION

Brand	Topform
User Height	163 cm to 178 cm
Seat Height	500 mm
Seat Depth	500 mm
Seat Width	480 mm
Seat To Top Of Back	640 mm
Chair Width	780 mm
Chair Height	1070 mm
Max SWL	155 kg
Weight capacity	100kg to 149kg