


# QUANTUM® 4FRONT<sup>2</sup>

## 4FRONT<sup>2</sup>

shown in *Military Green* with TRU-Balance<sup>®</sup> 4 seating system, Q-Logic 3 environmental controls, and Stealth<sup>®</sup> TRU-Comfort 2 cushion.



The 4Front<sup>®</sup> 2 is the next generation of the 4Front<sup>®</sup> Power Chair! The new Smart Traction Control (STC), along with Smooth Ride Suspension (SRS), provide stability for smooth obstacle transitions and unmatched driving performance. Additional new features include 4-pole motors, casted castor arms and a power-off switch.

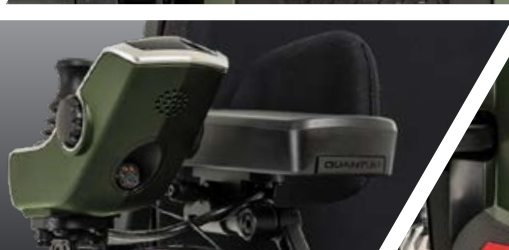

QUANTUM<sup>®</sup>  
#1 FOR REHAB POWER



quantumrehab.com

## Specifications

<b>DRIVE WHEELS</b>	Standard: 35.5 cm flat free Optional: 35.5 cm knobby/pneumatic
<b>CASTER WHEELS</b>	Rear: 22.8 cm flat free
<b>MAXIMUM SPEED</b> <sup>1</sup>	Up to 10 km/h Up to 5.6 km/h when elevated
<b>GROUND CLEARANCE</b>	15.2 cm at motor pivot
<b>TURNING RADIUS</b> <sup>2</sup>	65.7 cm
<b>OVERALL LENGTH</b> <sup>2</sup>	101.3 cm without front riggings
<b>BASE WIDTH</b>	61.4 cm
<b>SEATING SIZES</b>	
<b>TRU-Balance® 3 Static Seat</b>	<b>Width:</b> 12-22" <b>Depth:</b> 12-22"
<b>TRU-Balance® 3 Power Tilt</b>	<b>Width:</b> 12-22" <b>Depth:</b> 12-22"
<b>TRU-Balance® 3 Power Lift &amp; Tilt</b>	<b>Width:</b> 12-22" <b>Depth:</b> 12-22"
<b>SEAT-TO-FLOOR HEIGHTS</b>	
<b>TRU-Balance® 3 Static Seat</b> <sup>5</sup>	38.6-43.8 cm (not available with Synergy®)
<b>TRU-Balance® 3 Power Tilt</b> <sup>5</sup>	44.4-49.5 cm
<b>TRU-Balance® 3 Tilt &amp; iLevel®</b> <sup>5</sup>	45-50.1 cm
<b>BATTERY SIZE</b> <sup>3</sup>	Group 24 or Group 34
<b>BATTERY WEIGHT</b> <sup>7</sup>	Group 24: 23.3 kg Group 34: 20.3 kg
<b>AVAILABLE ELECTRONICS</b> <sup>8</sup>	120A Q-Logic 3e (standard) 120A Q-Logic 3 (optional)
<b>BATTERY CHARGER</b>	8A Off-board
<b>MOTOR PACKAGES</b>	4-pole, front-wheel drive
<b>WEIGHT CAPACITY</b> <sup>4</sup>	136 kg
<b>BASE WEIGHT</b>	58.9 kg

1) Speed varies with user weight, terrain type, battery charge, battery condition and tyre pressure.

2) Without front riggings

3) Battery capacity may change when adding power positioning systems.

4) Weight capacity is dictated by the selection of seating.

5) All measurements +/- 0.3 cm with 18" x 18" seat, fixed towers, 136 kg weight capacity and solid tyres (add 0.9 cm with pneumatic tyres). Widths may increase seat-to-floor height. Contact Customer Service for details. Listed measurements are in the lowest position; does not include elevating system unless specified.

6) Battery weight may vary +/- 1.3 kg, based on manufacturer.

7) All Quantum Rehab® units can be ordered with specialty controls.

8) Call for availability

9) Seat-to-floor heights are measured with the unit unloaded. Due to the 4Front's suspension system, seat-to-floor heights are approximately 0.6 cm lower when the unit is loaded.

## Standard Features

- 4-pole motors
- **NEW** Smart Traction Control (STC)
- Superior, automotive-grade SRS (Smooth Ride Suspension)
- Unmatched driving performance and ride quality
- Advanced stability for smooth obstacle transitions
- Casted castor arms
- Q-Logic 3e controller
- Increased speeds up to 10 km/h
- Power-off switch
- LED front and rear fender lights
- Quiet ride technology


## Options

- 25.4 cm of power adjustable safe seat elevation at speeds of up to 5.6 km/h
- **NEW** backup camera
- Trilogy vent tray <sup>8</sup>
- Compatible with TRU-Balance® 4 Seating System
- Compact vent tray <sup>8</sup>

## Colours



Red  
Corvette

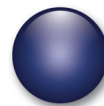

Midnight  
Blue



Electric  
Blue



Deep  
Purple



Lime  
Light



Back  
in Black



Military  
Green