

SUPINE TO PRONE IN ONE MOTION







Zing MPS, Supine, MPS TT(Tilt Table), Supine TT(Tilt Table), Vertical, & Prone

The Zing line of multi and single position standers were designed to fit the widest range of children and accommodate positioning challenges few others can. The Zing line features independent leg abduction up to 30 degrees on each model. The units can be configured simply, with padded strapping system to ease transfers and aid multi-user environments. But when required, individual positioning and support components can be added.



The Zing

True Multi-Position Stander

- The MPS base allows three positions without changing pads, flipping foot plates or moving the tray
- Provides a smooth transition from a flat-to-load transfer to supine, vertical and 20 degrees prone
- Frame will only go supine and upright with Prone Position Lockout engaged

True Hip Pivot Leg Abduction

- Pivot points are in line with the hips, promoting proper alignment in the hips, knees and feet
- · Feet maintain the angle with the rest of the legs no wedges needed
- · No need for abduction pad/pommel due to the innovative design
- Up to 30 degrees of abduction per leg
- · Promotes musculoskeletal development for children during critical growth years

Plenty of Growth

- Size 1: Fits children from infancy to 44", and up to 70lbs
- Size 2: Fits children from 40-60", and up to 154lbs
- · Can easily grow with the child and adapt to their changing needs
- · Many supports, pads and straps are available in multiple sizes

Remarkably Modular

- Can be configured as a simple stander or fit a child with the most complex needs
- Designed to grow with the child options can be added as needed
- Great for a multi-user facility or school

Advantage

Multi-Use Tray

- The MPS, Supine, MPS TT and Supine TT feature a multi-use tray that is angle adjustable and swings away either direction
- Tray can remain attached during transfers
- More activities are possible by swinging tray away during standing therapy

Unique Lift Mechanisms

- Dual controls are available allowing angle adjustment by foot pedal or hand control
- Gas Spring Lockout prevents an inadvertent change of position, especially important when around other children
- The Pow'r Up Lift option on Size 2 brings the user to the standing position with the touch of a button

Accommodating Contractures

- No special modifications are needed to accommodate children with contractures
- Using the multi-adjustable foot plates, multi-adjustable knee pads, and the independent calf pads, it is possible to adjust for contractures in one or both sides

Easy Adjustment

- Tool-free adjustment of all components
- · Every adjustment can be made with a simple turn of a knob or ratchet handle





Zing MPS, Supine, MPS TT(Tilt Table), Supine TT(Tilt Table), Vertical, & Prone

SUPINE TO PRONE

MPS Models Size 1 and Size 2










Health Benefits

- Improve/maintain range of motion
- Decrease joint/muscle contractures
- · Management of atrophy in the trunk and leg muscles
- Improve strength to trunk and lower extremities
- Decrease muscle spasms
- Promote proper muscle skeleton development
- Improve/maintain bone integrity/skeletal development
- · Lessen/manage the progression of scoliosis
- Manage pressure (ulcers) through changing positions
- Improve bowel function and regularity
- Aid in kidney and bladder functions
- Strengthen cardiovascular system and build endurance
- Improve circulation
- Reduce swelling in lower extremities

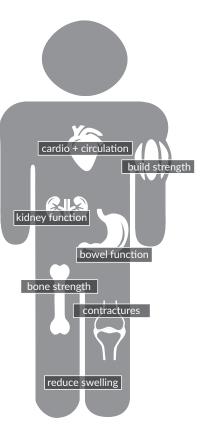

Zing Prone-Size 1

Hygienic Covers

- Available in 5 fun patterns for the pelvic, upper body, lateral, head and tray areas
- Prevents stander from getting soiled, especially during feeding
- Easily removed and washable



RAINBOW UNICORN



SEA CREATURES



FANCY PANTS ELEPHANTS



SPACE RACE



ALL SMILES ANIMALS

Frame Colors















Zing-Size 1 MPS, Supine, MPS TT(Tilt Table), Supine TT(Tilt Table), Vertical, & Prone




additional specifications at easystand.com

1 S J.com	MPS (Multi-Position Stander)	Supine	MPS TT (Tilt Table)
	The Zing MPS is the only multi-position stander that can seamlessly transition from supine to prone without having to transfer the child out and reposition foot plates, pads or trays.	The Zing Supine features a comfortable 32" transfer height and provides all the positioning and support features of the MPS, minus prone.	The Zing MPS TT is identical to the MPS from a positioning standpoint, but features larger, wider support pads and options for padded strapping, that can simplify transferring and securing children.
Weight Limit	70 lbs. (32kg)	70 lbs. (32kg)	70 lbs. (32kg)
Height Range	up to 44" (112cm)	up to 44" (112cm)	up to 44" (112cm)
Transfer Height	32"(81cm)	32" (81cm)	32" (81cm)
Pelvic Support to Foot Plate Range	7"-22" (18-56cm)	7"-22" (18-56cm)	7"-22" (18-56cm)
Upper Body Support Height Range	4" (10cm)	4" (10cm)	4" (10cm)
Knee Pad Depth Range	5" (13cm)	5" (13cm)	5" (13cm)
Weight of Basic Unit	55 lbs (25kg)	55 lbs (25kg)	55 lbs (25kg)
Frame Footprint	20"x28" (51x71cm)	20"x28" (51x71cm)	20"x28" (51x71cm)
30° Independent Leg Abduction	Available	Available	Available
Supine	0-90°	0-90°	0-90°
Prone	0°-20°	n/a	0°-20°
HCPCS	E0641	E0638	E0641
Positions	Supine, Vertical, Prone	Supine, Vertical	Supine, Vertical, Prone
Tray	Adjustable, Swing-Away	Adjustable, Swing-Away	Adjustable, Swing-Away
Support	Individual Positioning Components, Padded Strapping	Individual Positioning Components, Padded Strapping	Individual Positioning Components, Padded Strapping
Dual Control	Available	Available	Available
Gas Spring Lockout	Available	Available	Available
Transfer Aid	Flat-to-Load	Flat-to-Load	Flat-to-Load
Tool Free Adjustment	Yes	Yes	Yes

Supine TT (Tilt Table)	Vertical	Prone
The Zing Supine TT features the larger pads and padded strapping options of the MPS TT and positioning from flat to load supine to vertical.	The Zing Vertical is a single position stander that features a Padded Quick Strap to ease upright transfers. It's also one of the only vertical standers to offer anatomically correct independent leg abduction.	The Zing Prone is a single position stander but is one of the few to offer a flat-to-load transfer position at a comfortable height of 32". Prone positioning from 0-90 degrees.
70 lbs. (32kg)	70 lbs. (32kg)	70 lbs. (32kg)
up to 44" (112cm)	up to 44" (112cm)	up to 44" (112cm)
32" (81cm)	n/a	32" (81cm)
7"-22" (18-56cm)	7"-22" (18-56cm)	7"-22" (18-56cm)
4" (10cm)	4" (10cm)	4" (10cm)
5" (13cm)	5" (13cm)	5" (13cm)
55 lbs (25kg)	55 lbs (25kg)	55 lbs (25kg)
20"x28" (51x71cm)	20"x28" (51x71cm)	20"x28" (51x71cm)
Available	Available	Available
0-90°	n/a	n/a
n/a	n/a	0°-90°
E0638	E0638	E0638
Supine, Vertical	Vertical	Vertical, Prone
Adjustable, Swing-Away	Adjustable	Adjustable
Individual Positioning Components, Padded Strapping	Individual Positioning Components, Padded Strapping	Individual Positioning Components, Padded Strapping
Available	n/a	Available
Available	n/a	Available
Flat-to-Load	n/a	Flat-to-Load
Yes	Yes	Yes



FITTING THE ZING

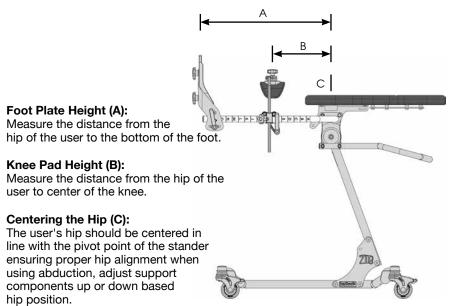

Zing MPS (Multi-Position Stander) Zing Supine Size 1-Options

The MPS (multi-position stander) base includes:

- Four Swivel Locking Casters
- Foot Operated Gas Spring
- Mast Receiver with Prone Position Lockout



PA5520 Zing MPS Base PA5522 Zing Supine Base

Comes with four swivel locking casters. foot operated gas spring, mast receiver. MPS base comes with prone position lockout.



PA5530 Mast

Mast is the base for positioning options including foot plates, knee supports, pelvic supports, upper body supports and head supports.



PA5532 Mast with Leg Abduction

Mast is the base for support and positioning components and allows for 0-30° of independent leg abduction. Multi-adjustable knee pads and multi-adjustable foot plates required for leg abduction.

PA5612 Dual Control

Allows both foot and hand control of the gas spring lift. Push handle PA5594 required.



PA5614 Gas Spring Lift Lockout Locks gas spring in any position from lying

down to fully upright, supine and prone positions.



PA5608 Directional Locking Caster Great for transporting stander. Directionally

locks one caster while leaving the other three free to swivel.



PA5536 Platform Foot Plate Adjustable in height. Foot holders and foot

straps can be added for more precise positioning of the feet.

The Supine base includes:

- Four Swivel Locking Casters
- Foot Operated Gas Spring
- Mast Receiver



Multi-Adjustable Foot Plates Foot plates are adjustable forward/aft +/-.75" (2cm) +extra 1.5" (3.8cm), plantar/dorsi +/- 20°, toe-in/toe-out, swivel 30° left/right, and in height. Foot holder sizes: PA5624-6"Lx3"W (15x8cm), PA5626-7.75"Lx3.25"W (20x8cm), PA5628-9.75"Lx4"W (25x10cm),

PA5540, PA5542, PA5544 Foot Holders for Platform Foot Plate

Use with the platform foot plate, gives added positioning of users feet. Foot holder sizes: PA5540-6"Lx3"W (15x8cm), PA5542-7.75"Lx3.25"W (20x8cm), PA5544-9.75"Lx4"W (25x10cm).

PA5546 Secure Foot Straps-8"L PT50082 Secure Foot Straps-10"L

Velcro® and D-ring adjustment holds feet securely in place at ankles and toes with multiple attachment slots for proper foot positioning. Straps measured in length over top of foot, with approximately 1.5" (4cm) of Velcro® overlap.

PA5548 Foot Straps-8"L PT50080 Foot Straps-10"L

Velcro® and D-ring adjustment holds feet in alignment with multiple attachment slots for proper foot positioning. Straps measured in length over top of foot, with approximately 1.5" (4cm) of Velcro® overlap.

PA5550 Calf/Knee Pad with Strap

PA5556, PA5558, PA5729, PA5730

4"x10" (10x25cm) pad attached to a 2" (5cm) Velcro® belt with a D-Ring closure, provides calf/knee stability for user. Fits circumference up to 24" (61cm). Gives support from the back of the calf/knee. Adjustable in height and fore/aft.



depth range of 5" (13cm).



Multi-Adjustable Knee Pads Multi-adjustable knee pads are adjustable in swivel, height +/-.75" (2cm), PA5556 inside width 2.5" (6cm), PA5558 inside width 3.25" (8cm), PA5729 inside width 4.25" (11cm) and PA5730 inside width 5" (13cm). With a

PA5560 Multi-Adjustable Calf Pads 3.5" PA5678 Multi-Adjustable Calf Pads 4.5"

Great for adapting to knee contractures. Gives support from the back of the calf. Adjustable in height with knee pad and fore/ aft 4" (10cm) independently. Must order multi-adjustable knee pads.





PA5562 Planar Pad-5"Hx7"W PA5564 Planar Pad-5"Hx9"W

Size 1

Cushioned pad covered with Neoprene®. Pad sizes: PA5562-5"Hx7"W (13x18cm), PA5564-5"Hx9"W (13x23cm).

70 lbs

(32 kg)

PA5566 Form to Fit Pad-5"Hx10"W PA5570 Form to Fit Pad-5"Hx12"W

Form to fit pad can be custom fit to the user. Pad sizes: PA5566-5"Hx10"W (13x25cm) PA5570-5"Hx12"W (13x30cm). Not recommended with leg abduction mast.

PA5574 Hip Supports

Keeps the hips in symmetrical alignment. Adjustable from 6"-10" (15-25cm) with fore/aft 4" (10cm). Pad size is 3.5"x3.5" (9x9cm). For planar pads only.

PA3011 Pelvic Positioning Belt

2" (5cm) wide belt has a Velcro® belt with a D-Ring closure, provides hip stability for user. Fits circumference up to 24" (61cm).

PA5610 Pelvic Padded Positioning Strap

5"x8" (13x20cm) pad attached to a 2" (5cm) belt, Velcro® belt with a D-Ring closure. Fits circumference up to 24" (61cm).

PA3028, PA3029 Hygienic Covers

Washable slip on cover available in five different patterns. For planar pads PA5562 and PA5564.

PA3030, PA3032 Hygienic Covers

Washable slip on cover available in five different patterns. For form to fit pads PA5566 and PA5570.











PA5562 Planar Pad-5"Hx7"W PA5564 Planar Pad-5"Hx9"W Cushioned pad covered with Neoprene®. Pad sizes: PA5562-5"Hx7"W (13x18cm), PA5564-5"Hx9"W (13x23cm).



PA5616 Planar Pad-7"Hx7"W PA5618 Planar Pad-7"Hx9"W Cushioned pad covered with Neoprene®. Pad sizes: PA5616-7"Hx7"W (18x18cm), PA5618-7"Hx9"W (18x23cm).



PA5568 Form to Fit Pad-5"Hx10"W PA5620 Form to Fit Pad-7"Hx10"W

Form to fit pad with built in lateral support that can be custom fit to the user. Pad size for PA5568 is 5"Hx10"W (13x25cm) and PA5620 is 7"Hx10"W (18x25cm).



PA5572 Form to Fit Pad-5"Hx12"W PA5622 Form to Fit Pad-7"Hx12"W

Form to fit pad with built in lateral support that can be custom fit to the user. Pad size for PA5572 is 5"Hx12"W (13x30cm) and PA5622 is 7"Hx12"W (18x30cm).



PA5576 Lateral Supports Keeps the upper torso in symmetrical alignment. Adjustable from 7-11"

(18-28cm) with a 3.5" (9cm) height range. Pad size is 3.5"Hx3.5"W (9x9cm). For planar pads only.



PA5578 Lateral Supports with Elbow Stop Works as a lateral support and elbow rest. Great for the supine position. For planar pads only.



PA5580 Lateral Supports w/Elbow Stop/Arm Rest Works as a lateral support, elbow rest and armrest. Great for the supine, upright and prone position. For planar pads only. Height and width adjustable.



PA5630 Elbow Stop with Arm Rest Great for loading and bringing to upright. Works with both planar and form to fit upper body supports. For planar and form to fit pads. Height and width adjustable.



PT50086 X-Style Chest Vest-9.5"Lx9"W PT50088 X-Style Chest Vest-11"Lx9.5"W Neoprene® with camlock adjustment provides anterior trunk support for child.





PA3028, PA3029, PA3045, PA3046 **Hygienic Covers** Washable slip on cover available in five different patterns. For planar pads PA5562, PA5564, PA5616, PA5618,

PA3031, PA3033, PA3047, PA3048 **Hygienic Covers** Washable slip on cover available in five

different patterns. For form to fit pads PA5568, PA5572, PA5620, PA5622.



PY5626 Head Support-5"Hx8"W PY5628 Head Support-6"Hx10"W Cushion size is 5"Hx8"W (13x20cm),

6"Hx10"W (15x25cm). Depth range is +1" to -2.75" (+2 to -7cm) from backrest. Maximum height from seat is 24" (61cm).

PA5586 Head Support-Form to Fit

Cushion size 5"Hx10"W (13x25cm), Depth range is -2.5" (6cm) to +1.5" (4cm) from backrest. Maximum height from top of upper body pad is 8" (20cm).

PA3035, PA3036 Hygienic Covers Washable slip on cover available in five different patterns. For PY5626 and PY5628.

PA3037 Hygienic Cover Washable slip on cover available in five different patterns. For PA5586.

PA5596 Black Molded Tray-11"Lx19"W

Tray is accessible in all positions and provides anterior support as user moves from supine to upright and prone positions. Adjustable in depth and height. Tray size is 11"Lx19"W (28x48cm).





Swing-away tray is accessible in all positions and provides anterior support as user moves from supine to upright and prone positions. Adjustable in depth, height and tray angle. Tray size is 11"Lx19"W (28x48cm).

PA5600 Black Molded Swing-Away Tray-11"Lx19"W

PA5602 Clear Swing Away Tray-11"Lx19"W

Swing-away tray is accessible in all positions and provides anterior support as user moves from supine to upright and prone positions. Adjustable in depth, height and trav angle. Trav size is 11"Lx19"W (28x48cm).

PA5604 Forearm Wings-Black Molded Tray

Forearm wings give additional area to rest arms when in standing position.



PA5606 Forearm Wings-Clear Tray Forearm wings give additional area to rest arms when in standing position.



PA3038 Hygienic Padded Cover

Washable slip on cover is padded and available in five different patterns.

PNG50042 Hand Grips

Provides additional arm and hand support and alignment. Easily movable for precise arm positioning. Attaches with suction cups. Clear tray required.

PA5594 Push Handle

Makes pushing the stander easier for the caregiver and also aids in lifting and lowering the child. Required for dual control.

P80809 Tool Pouch

Holds the assembly tools provided with the stander. Attaches to unit with Velcro® straps.



P80246 Angle Locator


Magnetic angle locater can be utilized for standing progress and leg abduction.









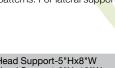





Zing MPS TT (Tilt Table) Zing Supine TT (Tilt Table) Size 1-Options

The MPS (multi-position stander) TT base includes:

- Four Swivel Locking Casters
- Foot Operated Gas Spring
- Mast Receiver with Prone Position Lockout



PA5522 Zing Supine TT Base PA5520 Zing MPS TT Base

Comes with four swivel locking casters. foot operated gas spring, mast receiver. MPS base comes with prone position lockout.



PA5632 Mast Mast is the base for support and positioning components.

PA5634 Mast with Leg Abduction Mast is the base for support and positioning components and allows for 0-30° of independent leg abduction. Multi-adjustable knee pads and multi-adjustable foot plates required for leg abduction.

PA5612 Dual Control

Allows both foot and hand control of the gas spring lift. Push handle PA5594 required.

PA5614 Gas Spring Lift Lockout

Locks gas spring in any position from lying down to fully upright, supine and prone positions.



PA5608 Directional Locking Caster

Great for transporting stander. Directionally locks one caster while leaving the other three free to swivel.

PA5536 Platform Foot Plate



Adjustable in height. Foot holders and foot straps can be added for more precise positioning of the feet.

The Supine TT base includes:

- Four Swivel Locking Casters
- Foot Operated Gas Spring
- Mast Receiver



Multi-Adjustable Foot Plates Foot plates are adjustable forward/aft +/-.75" (2cm) +extra 1.5" (3.8cm), plantar/dorsi +/- 20°, toe-in/toe-out, swivel 30° left/right, and in height. Foot holder sizes: PA5624-6"Lx3"W (15x8cm), PA5626-7.75"Lx3.25"W (20x8cm), PA5628-9.75"Lx4"W (25x10cm).

PA5540, PA5542, PA5544 Foot Holders for Platform Foot Plate

Use with the platform foot plate, gives added positioning of users feet. Foot holder sizes: PA5540-6"Lx3"W (15x8cm), PA5542-7.75"Lx3.25"W (20x8cm). PA5544-9.75"Lx4"W (25x10cm).

PA5546 Secure Foot Straps-8"L PT50082 Secure Foot Straps-10"L

Velcro® and D-ring adjustment holds feet securely in place at ankles and toes with multiple attachment slots for proper foot positioning. Straps measured in length over top of foot, with approximately 1.5" (4cm) of Velcro® overlap.

PA5548 Foot Straps-8"L PT50080 Foot Straps-10"L

PA5648-8"Hx6"W (20x15cm)

PA5644, PA5646, PA5648.

Velcro® and D-ring adjustment holds feet in alignment with multiple attachment slots for proper foot positioning. Straps measured in length over top of foot, with approximately 1.5" (4cm) of Velcro® overlap.





PA5550 Calf/Knee Pad with Strap



PA5660, PA5662, PA5664 Hygienic Covers

Washable slip on cover available in five

different patterns. For planar pads



PA5666, PA5668, PA5670 Knee Strap

Size 1

4"x6" (10x15cm) pad attached to a 1.5" (4cm) Velcro® belt with a D-Ring closure. Fits circumference up to 15" (38cm). Adjustable in height. For leg supports PA5644, PA5646, PA5648,

70 lbs

(32 kg)

0

PA5556, PA5558, PA5729, PA5730 Multi-Adjustable Knee Pads

Multi-adjustable knee pads are adjustable in swivel, height +/-.75" (2cm), PA5556 inside width 2.5" (6cm), PA5558 inside width 3.25" (8cm), PA5729 inside width 4.25" (11cm) and PA5730 inside width 5" (13cm). With a depth range of 5" (13cm).

PA5560 Multi-Adjustable Calf Pads 3.5" PA5678 Multi-Adjustable Calf Pads 4.5"

Great for adapting to knee contractures. Gives support from the back of the calf. Adjustable in height with knee pad and fore/aft -1"-4" (-2.5-10cm) independently. Must order multi-adjustable knee pads.


PA5638, PA5640, PA5642 Planar Pad

Cushioned pad covered with Neoprene®. PA5638-14"Hx14"W (36x36cm) PA5640-16"Hx14"W (41x36cm) PA5642-18"Hx14"W (46x36cm)

PA3063, PA3064, PA3065 Hygienic Covers

Washable slip on cover available in five different patterns. For planar pads PA5638, PA5640, PA5642.

PA5636 Hip Supports

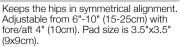

PA3062 Pelvic Positioning Strap

2" (5cm) wide belt has a Velcro® belt with a D-Ring closure, provides hip stability for user. Fits circumference up to 24" (61cm).









Pad sizes:







PA5672 Pelvic Padded Positioning Strap 5"x8" (13x20cm) pad attached to a 2" (5cm) belt. Fits circumference up to 24" (61cm).



PA5652 Lateral Supports-4"-8" Range PA5654 Lateral Supports-7"-11" Range Keeps the upper torso in symmetrical alignment. PA5652 Adjustable from 4"-8" (10-20cm). PA5654 Adjustable from 7"-11" (18-28cm) with a 3.5" (9cm) height range. Pad size is 3.5"Hx3.5"W (9x9cm).



PA5674 Upper Body Padded Positioning Strap 5"x8" (13x20cm) pad attached to a 2" (5cm) Velcro® belt with a D-Ring closure. Fits circumference up to 24" (61cm).



PT50086 X-Style Chest Vest-9.5"Lx9"W PT50088 X-Style Chest Vest-11"Lx9.5"W Neoprene® with camlock adjustment provides anterior trunk support for child.

PA3039 Hygienic Covers for Lateral Supports Washable slip on cover available in five

different patterns. For lateral supports

PA5652 and PA5654.



Swing-away tray is accessible in all positions and provides anterior support as user moves from supine to upright and prone positions. Adjustable in depth, height and tray angle. Tray size is 11"Lx19"W (28x48cm).

PA3066, PA3035, PA3036, PA3037

Washable slip on cover available in five

different patterns. For head supports

PA5650, PY5626, PY5628, PA5586,

PA5596 Black Molded Tray-11"Lx19"W

provides anterior support as user moves

positions. Adjustable in depth and height.

PA5600 Black Molded Swing-Away Tray-11"Lx19"W

Tray is accessible in all positions and

from supine to upright and prone

Tray size is 11"Lx19"W (28x48cm).

Swing-away tray is accessible in all

and prone positions. Adjustable in

11"Lx19"W (28x48cm).

positions and provides anterior support

depth, height and trav angle. Trav size is

PA5602 Clear Swing Away Tray-11"Lx19"W

as user moves from supine to upright

Hygienic Covers

PA5604 Forearm Wings-Black Molded Tray

arms when in standing position.



PA5650 Head Support-Planar Pad Cushioned Neoprene® pad for head support. Pad size 8"Hx8.5"W (21x22cm).




PY5626 Head Support-5"Hx8"W PY5628 Head Support-6"Hx10"W Cushion size is 5"Hx8"W (13x20cm), 6"Hx10"W (15x25cm). Depth range is +1" to -2.75" (+2 to -7cm) from backrest. Maximum height from seat is 24" (61cm).




PA5586 Head Support-Form to Fit Cushion size 5"Hx10"W (13x25cm), Depth range is -2.5" (6cm) to +1.5" (4cm) from backrest. Maximum height from top of upper body pad is 8" (20cm).



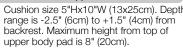

PA5676 Headrest Mounting Bracket For use of alternative manufacturer's head supports.



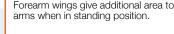

PA3038 Hygienic Padded Tray Cover

Washable slip on cover for tray is padded and available in five different patterns.

PNG50042 Hand Grips



Provides additional arm and hand support and alignment. Easily movable for precise arm positioning. Attaches with suction cups. Clear tray required.

PA5594 Push Handle

Makes pushing the stander easier for the caregiver and also aids in lifting and lowering the child. Required for dual control.





Holds the assembly tools provided with the stander. Attaches to unit with Velcro®

P80246 Angle Locator

Magnetic angle locater can be utilized for standing progress and leg abduction.

Forearm wings give additional area to rest





PA5606 Forearm Wings-Clear Tray Forearm wings give additional area to rest

Zing Vertical

Size 1-Options

The Vertical base includes:

- Four Swivel Locking Casters
- Mast Receiver



PA5524 Zing Vertical Base Comes with four swivel locking casters and a mast receiver.



PA5530 Mast Mast is the base for support and positioning components.



PA5532 Mast with Leg Abduction Mast is the base for support and positioning components and allows for 0-30° of independent leg abduction. Multi-adjustable knee pads and multi-adjustable foot plates required for leg abduction.



PA5608 Directional Locking Caster Great for transporting stander. Directionally locks one caster while leaving the other three free to swivel.




PA5680 Platform Foot Plate Adjustable in height. Foot holders and foot straps can be added for more precise positioning of the feet.



PA5624, PA5626, PA5628 Multi-Adjustable Foot Plates

Foot plates are adjustable forward/aft +/-.75" (2cm) +extra 1.5" (3.8cm), plantar/dorsi +/- 20°, toe-in/toe-out, swivel 30° left/right, and in height. Foot holder sizes: PA5624-6"Lx3"W (15x8cm), PA5626-7.75"Lx3.25"W (20x8cm), PA5628-9.75"Lx4"W (25x10cm).

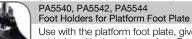

Use with the platform foot plate, gives added positioning of users feet. Foot holder sizes: PA5540-6"Lx3"W (15x8cm), PA5542-7.75"Lx3.25"W (20x8cm), PA5544-9.75"Lx4"W (25x10cm).











A M



4"x6" (10x15cm) pad attached to a 1.5" (4cm) Velcro® belt with a D-Ring closure. Fits circumference up to 15" (38cm). Adjustable in height. For PA5688, PA5690, PA5727 and PA5728.

PA5546 Secure Foot Straps-8"L

(4cm) of Velcro® overlap.

PA5548 Foot Straps-8"L PT50080 Foot Straps-10"L

overlap.

PT50082 Secure Foot Straps-10"L

for proper foot positioning. Straps

PA5550 Calf/Knee Pad with Strap

Adjustable in height and fore/aft.

PA5688 PA5690, PA5727, PA5728

Multi-Adjustable Knee Pads

4"x10" (10x25cm) pad attached to a 2"

provides calf/knee stability for user. Fits

circumference up to 24" (61cm). Gives

support from the back of the calf/knee.

(5cm) Velcro® belt with a D-Ring closure,



PA5566 Form to Fit Pad-5"Hx10"W PA5570 Form to Fit Pad-5"Hx12"W

Form to fit pad can be custom fit to the user. Pad sizes: PA5566-5"Hx10"W (13x25cm) PA5570-5"Hx12"W (13x30cm).



PA5574 Hip Supports


Keeps the hips in symmetrical alignment. Adjustable from 6"-10" (15-25cm) with fore/aft 4" (10cm). Pad size is 3.5"x3.5" (9x9cm). For planar pads only.

Size 1

70 lbs

∎(32 kg)



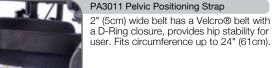

PA5610 Pelvic Padded Positioning Strap 5"x8" (13x20cm) pad attached to a 2" (5cm) Velcro® belt with a D-Ring closure. Fits circumference up to 24" (61cm).



PA3092 Quick Strap

2" (5cm) wide belt has a airline style buckle for quick closure, provides hip stability for user. Fits circumference up to 24" (61cm).


PA3093 Padded Quick Strap

5"x8" (13x20cm) pad attached to a 2" (5cm) belt has a airline style buckle for auick closure. Fits circumference up to 24" (61cm).



PA3028, PA3029, PA3045, PA3046 Hygienic Cover

Washable slip on cover available in five different patterns. For planar pads PA5562, PA5564, PA5616, PA5618.

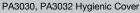

Washable slip on cover available in five different patterns. For form to fit pads PA5566 and PA5570.









PA5692, PA5562, PA5564, PA5616, PA5618 Planar Pad

Cushioned pad covered with Neoprene®. Pad sizes:

PA5692-5"Hx5"W (13x13cm) PA5562-5"Hx7"W (13x18cm) PA5564-5"Hx9"W (13x21cm) PA5616-7"Hx7"W (18x18cm) PA5618-7"Hx9"W (18x21cm)

PA5566 Form to Fit Pad-5"Hx10"W PA5570 Form to Fit Pad-5"Hx12"W

Form to fit pad with built in lateral support that can be custom fit to the user. Pad sizes are: PA5566-5"Hx10"W (13x25cm), PA5570-5"Hx12"W (13x30cm).



PA5694 Form to Fit Pad-5"Hx10"W PA5696 Form to Fit Pad-5"Hx12"W Form to fit pad with extra sternum positioning. Built in lateral support that can be custom fit to the user.. Pad sizes are: PA5694-5"Hx10"W (13x25cm), PA5696-5"Hx12"W (13x30cm).



PA5698 Form to Fit Pad-7"Hx10"W PA5700 Form to Fit Pad-7"Hx12"W Form to fit pad with extra sternum positioning. Built in lateral support that can be custom fit to the user. Pad sizes are: PA5698-7"Hx10"W (18x25cm), PA5700-7"Hx12"W (18x30cm).



PA5702 Lateral Supports-4"-8"Range PA5704 Lateral Supports-7"-11"Range Keeps the upper torso in symmetrical alignment. PA5702 Adjustable from 4"-8" (10-20cm), PA5704 Adjustable from 7"-11" (18-28cm) with a 3.5" (9cm) height range. Pad size is 3.5"Hx3.5"W (9x9cm). For planar pads only.




PA5610 Upper Body Padded Positioning Strap 5"x8" (13x20cm) pad attached to a 2" (5cm) Velcro® belt with a D-Ring closure. Fits circumference up to 24" (61cm).



PA3039 Hygienic Covers for Lateral Supports Washable slip on cover available in five different patterns. For lateral supports PA5702 and PA5704.



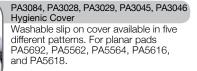

PA3030, PA3032 Hygienic Cover Washable slip on cover available in five different patterns. For form to fit pads PA5566 and PA5570.







PA5712 Black Molded Tray-11"Lx19"W Trav provides anterior support. Adjustable in depth and height. Tray size is 11"Lx19"W (28x48cm).

PA5710 Hygienic Cover

PA5706 and PA5708.

PA3085, PA3086, PA3087, PA3088

PA5706 Face Aperture-5"H

PA5708 Face Aperture-7"H

Washable slip on cover available in five

Washable slip on cover available in five

different patterns. For face apertures

different patterns. For form to fit pads

Hygienic Cover

PA5716 Multi-Adjustable Black Molded Tray-11"Lx19"W

Tray provides anterior support. Adjustable in depth, height and tray angle. Tray size is 11"Lx19"W (28x48cm).


PA5718 Multi-Adjustable Clear Tray-11"Lx19"W

Tray provides anterior support. Adjustable in depth, height and tray angle. Tray size is 11"Lx19"W (28x48cm).

PA5604 Forearm Wings-Black Molded Tray

Forearm wings give additional area to rest arms when in standing position.

PA5606 Forearm Wings-Clear Trav

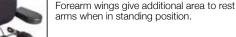

PA3038 Hygienic Padded Cover




Washable slip on cover for tray is padded and available in five different patterns.



PNG50042 Hand Grips

Provides additional arm and hand support and alignment. Easily movable for precise arm positioning. Attaches with suction cups. Clear tray required.

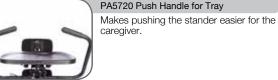

P80809 Tool Pouch

Holds the assembly tools provided with the stander. Attaches to unit with Velcro® straps.



P80246 Angle Locator

Magnetic angle locater can be utilized for leg abduction.







Zing Prone

Size 1-Options

The Prone base includes:

- Four Swivel Locking Casters
- Foot Operated Gas Spring
- Mast Receiver



PA5522 Zing Prone Base Comes with four swivel locking casters, foot operated gas spring, and a mast receiver.



PA5530 Mast Mast is the base for support and

positioning components.




PA5532 Mast with Leg Abduction Mast is the base for support and positioning components and allows for 0-30° of independent leg abduction.

Multi-adjustable knee pads and

multi-adjustable foot plates required for leg abduction.

PA5725 Dual Control

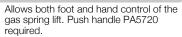



PA5614 Gas Spring Lift Lockout Locks gas spring in any position from lying

down to fully upright and supine positions.




PA5608 Directional Locking Caster

Great for transporting stander. Directionally locks one caster while leaving the other three free to swivel.




page 13

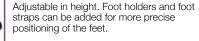

PA5680 Platform Foot Plate



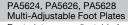

Foot plates are adjustable forward/aft +/-.75" (2cm) +extra 1.5" (3.8cm), plantar/dorsi +/- 20°, toe-in/toe-out, swivel 30° left/right, and in height. Foot holder sizes: PA5624-6"Lx3"W (15x8cm), PA5626-7.75"Lx3.25"W (20x8cm), PA5628-9.75"Lx4"W (25x10cm),

PA5540, PA5542, PA5544 Foot Holders for Platform Foot Plate

Use with the platform foot plate, gives added positioning of users feet. Foot holder sizes: PA5540-6"Lx3"W (15x8cm). PA5542-7.75"Lx3.25"W (20x8cm), PA5544-9.75"Lx4"W (25x10cm).

PA5546 Secure Foot Straps-8"L PT50082 Secure Foot Straps-10"L

Velcro® and D-ring adjustment holds feet securely in place at ankles and toes with multiple attachment slots for proper foot positioning. Straps measured in length over top of foot, with approximately 1.5" (4cm) of Velcro® overlap.

PA5548 Foot Straps-8"L PT50080 Foot Straps-10"L

Velcro® and D-ring adjustment holds feet in alignment with multiple attachment slots for proper foot positioning. Straps measured in length over top of foot, with approximately 1.5" (4cm) of Velcro® overlap.

PA5550 Calf/Knee Pad with Strap

4"x10" (10x25cm) pad attached to a 2" (5cm) Velcro® belt with a D-Ring closure, provides calf/knee stability for user. Fits circumference up to 24" (61cm). Gives support from the back of the calf/knee. Adjustable in height and fore/aft.

PA5688 PA5690, PA5727, PA5728 Multi-Adjustable Knee Pads

Multi-adjustable knee pads are adjustable in swivel, height +/-.75" (2cm), PA5688 inside width 2.5" (6cm), PA5690 inside width 3.25" (8cm), PA5727 inside width 4.25" (11cm) and PA5728 inside width 5" (13cm). With a depth range of 5" (13cm).

PA5726 Knee Strap

4"x6" (10x15cm) pad attached to a 1.5" (4cm) Velcro® belt with a D-Ring closure. Fits circumference up to 15" (38cm). Adjustable in height. For PA5688, PA5690, PA5727 and PA5728.



PA5566 Form to Fit Pad-5"Hx10"W PA5570 Form to Fit Pad-5"Hx12"W

Form to fit pad can be custom fit to the user. Pad sizes: PA5566-5"Hx10"W (13x25cm)PA5570-5"Hx12"W (13x30cm).

Cushioned pad covered with Neoprene®.

Pad sizes: PA5562-5"Hx7"W (13x18cm),

Size 1

PA5562 Planar Pad-5"Hx7"W

PA5564 Planar Pad-5"Hx9"W

PA5564-5"Hx9"W (13x23cm).

70 lbs

(32 kg)


PA5574 Hip Supports

Keeps the hips in symmetrical alignment. Adjustable from 6"-10" (15-25cm) with fore/aft 4" (10cm). Pad size is 3.5"x3.5" (9x9cm). For planar pads only.

PA3011 Pelvic Positioning Strap

2" (5cm) wide belt has a Velcro® belt with a D-Ring closure, provides hip stability for user. Fits circumference up to 24" (61cm).

PA5610 Pelvic Padded Positioning Strap

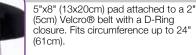



PA3092 Quick Strap

2" (5cm) wide belt has a airline style buckle for quick closure, provides hip stability for user. Fits circumference up to 24" (61cm).

PA3093 Padded Quick Strap

5"x8" (13x20cm) pad attached to a 2" (5cm) belt has a airline style buckle for quick closure. Fits circumference up to 24" (61cm).












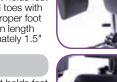



PA3028, PA3029, PA3045, PA3046 Hygienic Cover Washable slip on cover available in five different patterns. For planar pads PA5562, PA5564, PA5616, PA5618,



Sana P

PA3030, PA3032 Hygienic Cover

Washable slip on cover available in five different patterns. For form to fit pads PA5566 and PA5570.

PA5692, PA5562, PA5564, PA5616, PA5618

Cushioned pad covered with Neoprene®. Pad

positioning. Built in lateral support that can

be custom fit to the user. Pad sizes are:

PA5566 5"Hx10"W (13x25cm), PA5570

PA5692-5"Hx5"W (13x13cm) PA5562-5"Hx7"W (13x13cm) PA5562-5"Hx7"W (13x18cm) PA5564-5"Hx9"W (13x21cm) PA5616-7"Hx7"W (18x18cm)

PA5618-7"Hx9"W (18x21cm) PA5566 Form to Fit Pad-5"Hx10"W

5"Hx12"W (13x30cm).


PA5570 Form to Fit Pad-5"Hx12"W

Form to fit pad with extra sternum

Planar Pad

sizes



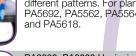

Hygienic Cover



PA3084, PA3028, PA3029, PA3045, PA3046

PA3030, PA3032 Hygienic Cover

Washable slip on cover available in five different patterns. For form to fit pads PA5566 and PA5570.

PA3085, PA3086, PA3087, PA3088 Hygienic Cover

Washable slip on cover available in five different patterns. For form to fit pads PA5694, PA5696, PA5698 and PA5700.

Pad size is 9"Hx3"W (23x8cm). With a 5"

(13cm) or 7" (18cm) range from the upper


body support, individual face aperture

pads are adjustable. Supports the head.

PA5706 Face Aperture-5"H

PA5708 Face Aperture-7"H



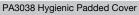

PA5606 Forearm Wings-Clear Tray

arms when in standing position.

Washable slip on cover for tray is padded and available in five different patterns.

Forearm wings give additional area to rest

PNG50042 Hand Grips

Provides additional arm and hand support and alignment. Easily movable for precise arm positioning. Attaches with suction cups. Clear tray required.



PA5720 Push Handle for Trav Makes pushing the stander easier for the caregiver and also aids in lifting and

lowering the child. Required for dual control. Not available with PA5712.



PA5694 Form to Fit Pad-5"Hx10"W PA5696 Form to Fit Pad-5"Hx12"W Form to fit pad with extra sternum

positioning. Built in lateral support that can be custom fit to the user. Pad sizes are: PA5694-5"Hx10"W (13x25cm), PA5696-5"Hx12"W (13x30cm).



PA5698 Form to Fit Pad-7"Hx10"W PA5700 Form to Fit Pad-7"Hx12"W Form to fit pad with extra sternum

positioning. Built in lateral support that can be custom fit to the user. Pad sizes are: PA5698-7"Hx10"W (18x25cm), PA5700-7"Hx12"W (18x30cm).



PA5702 Lateral Supports-4"-8"Range PA5704 Lateral Supports-7"-11"Range Keeps the upper torso in symmetrical alignment, PA5702 Adjustable from 4"-8" (10-20cm), PA5704 Adjustable from 7"-11" (18-28cm) with a 3.5" (9cm) height range. Pad size is 3.5"Hx3.5"W (9x9cm). For planar pads only.

PA5610 Upper Body Padded Positioning Strap 5"x8" (13x20cm) pad attached to a 2"

(5cm) Velcro® belt with a D-Ring closure. Fits circumference up to 24" (61cm).

PA3039 Hygienic Cover for Lateral Supports



Washable slip on cover available in five different patterns. For lateral supports PA5702 and PA5704.



PA5716 Multi-Adjustable Black Molded Tray-11"Lx19"W

is 11"Lx19"W (28x48cm).

Adjustable in depth and height. Tray size

Tray provides anterior support. Adjustable in depth, height and tray angle. Tray size is 11"Lx19"W (28x48cm).

PA5718 Multi-Adjustable Clear Tray-11"Lx19"W

Trav provides anterior support. Adjustable in depth, height and tray angle. Tray size is 11"Lx19"W (28x48cm).

PA5604 Forearm Wings-Black Molded Tray

Forearm wings give additional area to rest arms when in standing position.








P80246 Angle Locator


Magnetic angle locater can be utilized for standing progress and leg abduction.

P80809 Tool Pouch Holds the assembly tools provided with the stander. Attaches to unit with Velcro® straps.

page 14

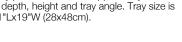



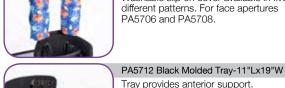

Zing-Size 2 MPS, Supine, MPS TT(Tilt Table), Supine TT(Tilt Table), & Prone



additional specifications at easystand.com

2 S d.com	MPS (Multi-Position Stander)	Supine	MPS TT (Tilt Table)
	The Zing MPS is the only multi-position stander that can seamlessly transition from supine to prone without having to transfer the child out and reposition foot plates, pads or trays.	The Zing Supine features a comfortable 31 ⁺ transfer height and provides all the positioning and support features of the MPS, minus prone.	The Zing MPS TT is identical to the MPS from a positioning standpoint, but features larger, wider support pads and options for padded strapping, that can simplify transferring and securing children.
Weight Limit	154 lbs. (70kg)	154 lbs. (70kg)	154 lbs. (70kg)
Height Range	40"-60" (102-152cm)	40"-60" (102-152cm)	40"-60" (102-152cm)
Transfer Height	31" (79cm)	31" (79cm)	31" (79cm)
Hip to Foot Plate Range	34" (86cm)	34" (86cm)	34" (86cm)
Upper Body Support Height Range	6" (15cm)	6" (15cm)	6" (15cm)
Knee Pad Depth Range	7" (18cm)	7" (18cm)	7" (18cm)
Weight of Basic Unit	82 lbs (37kg)	82 lbs (37kg)	82 lbs (37kg)
Frame Footprint	26.75"x39.75" (68x101cm)	26.75"x39.75" (68x101cm)	26.75"x39.75" (68x101cm)
30° Independent Leg Abduction	Available	Available	Available
Supine	0-90°	0-90°	0-90°
Prone	0°-20°	n/a	0°-20°
HCPCS	E0641	E0638	E0641
Positions	Supine, Vertical, Prone	Supine, Vertical	Supine, Vertical, Prone
Tray	Adjustable, Swing-Away	Adjustable, Swing-Away	Adjustable, Swing-Away
Support	Individual Positioning Components, Padded Strapping	Individual Positioning Components, Padded Strapping	Individual Positioning Components, Padded Strapping
Dual Control	Available	Available	Available
Gas Spring Lockout	Available	Available	Available
2 Transfer Aid	Flat-to-Load	Flat-to-Load	Flat-to-Load
Tool Free Adjustment	Yes	Yes	Yes
Pow'r Up Lift Option	Yes	Yes	Yes

Zing MPS- Size 2

page 15

Supine TT (Tilt Table)	Prone	
The Zing Supine TT features the larger pads and padded strapping options of the MPS TT and positioning from flat to load supine to vertical.	The Zing Prone is a single position stander but is one of the few to offer a flat-to-load transfer position at a comfortable height of 31". Prone positioning from 0-90 degrees.	
154 lbs. (70kg)	154 lbs. (70kg)	
40"-60" (102-152cm)	40"-60" (102-152cm)	
31" (79cm)	31" (79cm)	
34" (86cm)	34" (86cm)	
6" (15cm)	6" (15cm)	
7" (18cm)	7" (18cm)	
82 lbs (37kg)	82 lbs (37kg)	
26.75"x39.75" (68x101cm)	26.75"x39.75" (68x101cm)	
Available	Available	
0-90°	n/a	
n/a	0°-90°	
E0638	E0638	
Supine, Vertical	Vertical, Prone	
Adjustable, Swing-Away Individual Positioning Components, Padded Strapping	Adjustable Individual Positioning Components, Padded Strapping	
Available	Available	
Available	Available	
Flat-to-Load	Flat-to-Load	
Yes	Yes	
Yes	Yes	



FITTING THE ZING

Foot Plate Height (A):
Measure the distance from the hip of the user to the bottom of the foot.
Knee Pad Height (B):
Measure the distance from the hip of the user to center of the knee.
Centering the Hip (C):
The user's hip should be centered in line with the pivot point of the stander ensuring proper hip alignment when using abduction, adjust support components up or down based hip position.

Zing MPS (Multi-Position Stander) Zing Supine Size 2-Options

The MPS (multi-position stander) base includes:

- Four Swivel Locking Casters
- Foot Operated Gas Spring or Pow'r Up Lift
- Mast Receiver with Prone Position Lockout



PB5500 Zing MPS Base Gas Spring Lift PB5502 Zing Supine Base Gas Spring Lift Comes with four swivel locking casters. foot operated gas spring, mast receiver. MPS base comes with prone position lockout.



PB5504 Zing MPS Base Pow'r Up Lift PB5506 Zing Supine Base Pow'r Up Lift Comes with four swivel locking casters, Pow'r Up Lift, mast receiver. MPS base comes with prone position lockout.

PB5512 Mast



Mast is the base for positioning options including foot plates, knee supports, pelvic supports, upper body supports and head



PB5514 Mast with Leg Abduction

Mast is the base for support and positioning components and allows for 0-30° of independent leg abduction. Multi-adjustable knee pads and multi-adjustable foot plates required for leg abduction.

PB5516 Dual Control



Allows both foot and hand control of the gas spring lift. Push handle PB5560 required.



PA5614 Gas Spring Lift Lockout

Locks gas spring in any position from lying down to fully upright, supine and prone positions.



PB5510 Directional Locking Caster Great for transporting stander. Directionally locks one caster while leaving the other three free to swivel.

The Supine base includes:

- Four Swivel Locking Casters
- Foot Operated Gas Spring or Pow'r Up Lift

PB5518 Platform Foot Plate

positioning of the feet. Platform

16.75"Wx11.5"L (29-43cm).

Mast Receiver





PB5520, PB5522, PB5524 Multi-Adjustable Foot Plates Foot plates are adjustable forward/aft +/-.75" (2cm) +extra 1.5" (3.8cm), plantar/dorsi +/- 20°, toe-in/toe-out, swivel 30° left/right, and in height. Foot holder sizes: PB5520 7.75"Lx3.25"W (20x8cm), PB5522 9.75"Lx4"W (25x10cm). PB5524 11.75"Lx5"W (30x13cm).

Adjustable in height. Foot holders and foot

straps can be added for more precise

PA5542, PA5544, PB5526 Foot Holders for Platform Foot Plate Use with the platform foot plate, provides added positioning of users feet. Foot holder sizes: PA5542 7.75"Lx3.25"W (20x8cm),PA5544 9.75"Lx4"W (25x10cm), PB5526 11.75"Lx5"W (30x13cm).

PT50082 Secure Foot Straps-10"L PY5638 Secure Foot Straps-15"L PNG50506 Secure Foot Straps-21"L Velcro® and D-ring adjustment holds feet securely in place at ankles and toes with multiple attachment slots for proper foot positioning. Straps measured in length over top of foot, with approximately 1.5" (4cm) of Velcro® overlap.

PT50080 Foot Straps-10"L PY5636 Foot Straps-15"L PNG50505 Foot Straps-21"L

Velcro® and D-ring adjustment holds feet in alignment with multiple attachment slots for proper foot positioning. Straps measured in length over top of foot, with approximately 1.5" (4cm) of Velcro® overlap

PB5528 Calf/Knee Pad with Strap

17.5"x5.5" (44x14cm) pad attached to a 2" (5cm) Velcro® belt with a D-Ring closure, provides calf/knee stability for user. Fits circumference up to 28"-38" (71-97cm). Gives support from the back of the calf/ knee. Adjustable in height and fore/aft.

PB5530, PB5532, PB5534 Multi-Adjustable Knee Pads

Multi-adjustable knee pads are adjustable in swivel, height +/-.75" (2cm), PA5530 inside width 3.25" (8cm), PA5532 inside width 4.25" (11cm) and PB5534 inside width 5" (13cm). With a depth range of 7" (18cm).





3 1 4

PA5678 Multi-Adjustable Calf Pads 4.5" PB5536 Multi-Adjustable Calf Pads 5.5"

Size 2

154 lbs (70 kg)

Great for adapting to knee contractures. Gives support from the back of the calf. Adjustable in height with knee pad and fore/aft 4" (10cm) independently. Must order multi-adjustable knee pads.

PB3000 Planar Pad 10"Hx11"W

PB3002 Planar Pad 10"Hx13"W Cushioned pad covered with Neoprene®. Pad sizes: PB3000-10"Hx11"W (25x28cm)

PB3002-10"Hx13"W (25x33cm)

PB3004 Form to Fit Pad 10"Hx11"W PB3006 Form to Fit Pad 10"Hx13"W


Form to fit pad can be custom fit to the user. Pad sizes: PB3004 10"Hx11"W (25x28cm) PB3006 10"Hx13"W (25x33cm). Not recommended with lea abduction mast.

PB5538, PB5540 Hip Supports

Keeps the hips in symmetrical alignment. Adjustable from PB5538 7"-13" (18-33cm) and PB5540 12"-18" (30-46cm) with fore/aft 4" (10cm). Pad size is 4.5"x4" (11x10cm). For planar pads only.

PT30066 Pelvic Positioning Strap

2" (5cm) wide belt has a Velcro® belt with a D-Ring closure, provides hip stability for user. Fits circumference up to 35" (89cm).

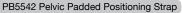

6"x10" (15x25cm) pad attached to a 2" (5cm) Velcro® belt with a D-Ring closure. Fits circumference up to 35" (89cm).

PB3025, PB3026 Hygienic Cover

Washable slip on cover available in five different patterns. For planar pads PB3000 and PB3002.











PB3027, PB3028 Hygienic Cover Washable slip on cover available in five different patterns. For form to fit pads PB3004 and PB3006.



PB3001 Planar Pad-8"Hx11"W PB3003 Planar Pad-8"Hx13"W Cushioned pad covered with Neoprene®. Pad size for PB3001 is 8"Hx11"W (20x28cm) and PB3003 is 8"Hx13"W (20x33cm).



PB3005 Form to Fit Pad-8"Hx11"W PB3007 Form to Fit Pad-8"Hx13"W Form to fit pad with built in lateral support that can be custom fit to the user. Pad size for PB3005 is 8"Hx11"W (20x28cm) and PB3007 is 8"Hx13"W (20x33cm).



PB5544, PB5546 Lateral Supports

Keeps the upper torso in symmetrical alignment. Adjustable from PB5544 7"-11"W with PB3001 and 9"-13"W with PB3003 and PB5546 10"-14"W with PB3001 and 12"-16"W with PB3003 with a 3.5" (9cm) height range. Pad size is 4"Hx6"W (10x15cm). For planar pads only.



PB5548 Lateral Supports with Elbow Stop

Works as a lateral support and elbow rest. Great for the supine position. Range: 10"-17" (25-43cm) with planar pad B3001 and 12"-19" (30-48cm) with planar pad B3003. For planar pads only.



PB5550 Lateral Supports w/Elbow Stop/Arm Rest Works as a lateral support, elbow rest and armrest. Great for the supine, upright and prone position. For planar pads only. Height and width adjustable.



PB5552 Elbow Stop with Arm Rest Great for loading and bringing to upright. Works with both planar and form to fit upper body supports. For planar pads only. Height and width adjustable.



PNG30170 X-Style Chest Vest-9.5"Lx9"W PNG30056 X-Style Chest Vest-11"Lx9.5"W Neoprene® with camlock adjustment provides anterior trunk support for child.




PNG30171 Y-Style Chest Vest-17.5"Lx11.5"W PNG30021 Y-Style Chest Vest-20.5"Lx12.5"W Neoprene® with camlock adjustment provides anterior trunk support for child.







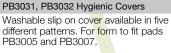

PB5542 Upper Body Padded Positioning Strap

6"x10" (15x25cm) pad attached to a 2"

closure. Fits circumference up to 35"

Washable slip on cover available in five

different patterns. For planar pads PB3301

(5cm) Velcro® belt with a D-Ring

PB3029, PB3030 Hygienic Cover

(89cm).

and PB3003.

PB3049 Hygienic Cover for PB5544/PB5546

Washable slip on cover available in five different patterns. For PB5544 and PB5546.

PY5626 Head Support 5"Hx8"W PY5628 Head Support 6"Hx10"W PY5630 Head Support 7"Hx14"W

Cushion sizes are: PY5626 5"Hx8"W (13x20cm). PY5628 6"Hx10"W (15x25cm), PY5630 7"Hx114"W (18x36cm). Depth range is +1" to -2.75" (+2 to -7cm) from backrest. Maximum height from top of upper body pad is 10" (25cm).

PA5586 Head Support-Form to Fit

Cushion size 5"Hx10"W (13x25cm). Depth range is -2.5" (6cm) to +1.5" (4cm) from backrest. Maximum height from top of upper body pad is 10" (25cm).

PA3035, PA3036, PB3033, PA3037 Hygienic Cover

Washable slip on cover available in five different patterns. For PY5626, PY5628 PY5630 and PA5586.

PB5554 Black Molded Tray-21"Lx24"W

Trav is accessible in all positions and provides anterior support as user moves from supine to upright and prone positions. Adjustable in depth and height. Tray size is 21"Lx24"W (53x61cm).

PB5556 Black Molded Swing-Away Tray-21"Lx24"W

Swing-away tray is accessible in all positions and provides anterior support as user moves from supine to upright and prone positions. Adjustable in depth, height and tray angle. Tray size is 21"Lx24"W (53x61cm).



PB5558 Clear Swing Away Tray-21"Lx24"W

Swing-away tray is accessible in all positions and provides anterior support as user moves from supine to upright and prone positions. Adjustable in depth, height and tray angle. Tray size is 21"Lx24"W (53x61cm).

PT50156 Elbow Pad-9.25" Cutout

Additional padded support surface on the black molded or clear tray. 9.25" (24cm) cutout width. Overall width 21.5" (55cm), pad extends back 3.75" (10cm) posterior from current table edge for additional upper body/extremity support.

PT50160 Elbow Pad-11" Cutout

Additional padded support surface on the black molded or clear tray. 11" (28cm) cutout width. Overall width 21.5" (55cm). pad extends back 3.75" (10cm) posterior from current table edge for additional upper body/extremity support.

PY5582 Elbow Pad-13" Cutout

Additional padded support surface on the black molded or clear tray. 13" (33cm) cutout width. Overall width 21.5" (55cm). pad extends back 3.75" (10cm) posterior from current table edge for additional upper body/extremity support.

PNG30295 Hygienic Padded Cover

Washable slip on cover is padded and available in five different patterns.

PNG50042 Hand Grips

Provides additional arm and hand support and alignment. Easily movable for precise arm positioning. Attaches with suction cups. Clear tray required.

PB5560 Push Handle

Makes pushing the stander easier for the caregiver and also aids in lifting and lowering the child. Required for dual control.

P80809 Tool Pouch

Holds the assembly tools provided with the stander. Attaches to unit with Velcro® straps.

P80246 Angle Locator

Magnetic angle locater can be utilized for standing progress and leg abduction.



















Zing MPS TT (Tilt Table) Zing Supine TT (Tilt Table) Size 2-Options

The MPS (multi-position stander) TT base includes:

- Four Swivel Locking Casters
- Foot Operated Gas Spring or Pow'r Up Lift
- Mast Receiver with Prone Position Lockout



PB5562 Zing MPS TT Base Gas Spring Lift PB5564 Zing Supine TT Base Gas Spring Lift Comes with four swivel locking casters. foot operated gas spring, mast receiver. MPS base comes with prone position lockout.



PB5566 Zing MPS TT Base Pow'r Up Lift PB5568 Zing Supine TT Base Pow'r Up Lift Comes with four swivel locking casters, Pow'r Up Lift, mast receiver. MPS base comes with prone position lockout.



PB5512 Mast


Mast is the base for positioning options including foot plates, knee supports, pelvic supports, upper body supports and head supports.



PB5514 Mast with Leg Abduction

Mast is the base for support and positioning components and allows for 0-30° of independent leg abduction. Multi-adjustable knee pads and multi-adjustable foot plates required for leg abduction.


PB5516 Dual Control

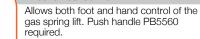



PA5614 Gas Spring Lift Lockout

Locks gas spring in any position from lying down to fully upright, supine and prone positions.

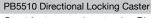

Great for transporting stander. Directionally locks one caster while leaving the other three free to swivel.

The Supine TT base includes:

- Four Swivel Locking Casters
- Foot Operated Gas Spring or Pow'r Up Lift
- Mast Receiver



PB5518 Platform Foot Plate Adjustable in height. Foot holders and foot straps can be added for more precise positioning of the feet.



(25x10cm). PB5524 11.75"Lx5"W (30x13cm). PA5542, PA5544, PB5526 Foot Holders for Platform Foot Plate

Use with the platform foot plate, gives added positioning of users feet. Foot holder sizes: PA5542 7.75"Lx3.25"W (20x8cm), PA5544 9.75"Lx4"W (25x10cm) PB5526 11.75"Lx5"W (30x13cm).

PT50082 Secure Foot Straps-10"L PY5638 Secure Foot Straps-15"L PNG50506 Secure Foot Straps-21"L Velcro® and D-ring adjustment holds feet securely in place at ankles and toes with multiple attachment slots for proper foot positioning. Straps measured in length over top of foot, with approximately 1.5" (4cm) of Velcro® overlap.

PT50080 Foot Straps-10"L PY5636 Foot Straps-15"L PNG50505 Foot Straps-21"L

Velcro® and D-ring adjustment holds feet in alignment with multiple attachment slots for proper foot positioning. Straps measured in length over top of foot, with approximately 1.5" (4cm) of Velcro®

Cushioned Neoprene® pad for leg PB5570- 9"Hx7.5"W (19x23cm) PB5572-12"Hx7.5"W (19x30cm) PB5574-15"Hx7.5"W (19x38cm)

PB5600, PB5602, PB5604 Hygienic Covers

Washable slip on cover available in five different patterns. For planar pads PB5570, PB5572, PB5574.



PB5528 Calf/Knee Pad with Strap

17.5"x5.5" (44x14cm) pad attached to a 2" (5cm) Velcro® belt with a D-Ring closure. provides calf/knee stability for user. Fits circumference up to 28"-38" (71-97cm). Gives support from the back of the calf/ knee. Adjustable in height and fore/aft.

Size 2

154 lbs (70 kg)

PB5576, PB5578, PB5580 Knee Strap

4"x6" (10x15cm) pad attached to a 1.5" (4cm) Velcro® belt with a D-Ring closure. Fits circumference up to 19" (48cm). Adjustable in height. For leg supports PB5570, PB5572, PB5574.

PB5530, PB5532, PB5534 Multi-Adjustable Knee Pads

Multi-adjustable knee pads are adjustable in swivel, height +/-.75" (2cm), PA5530 inside width 3.25" (8cm), PA5532 inside width 4.25" (11cm) and PB5534 inside width 5" (13cm). With a depth range of 7" (18cm).

PA5678 Multi-Adjustable Calf Pads 4.5" PB5536 Multi-Adjustable Calf Pads 5.5" Great for adapting to knee contractures. Gives support from the back of the calf.

Adjustable in height with knee pad and fore/aft 4" (10cm) independently. Must order multi-adjustable knee pads.

PB3012, PB3015, PB3018 Planar Pad


Cushioned pad covered with Neoprene®. Pad sizes:

PB3012-19"Hx16"W (48x41cm) PB3015-22"Hx16"W (56x41cm) PB3018-25"Hx16"W (64x51cm)



PB3035, PB3036, PB3037 Hygienic Cover

Washable slip on cover available in five different patterns. For planar pads PB3012, PB3015, PB3018.

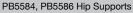

Keeps the hips in symmetrical alignment. Adjustable from PB5538 7"-13" (18-33cm) and PB5540 10"-16" (25-41cm) with fore/aft 4" (10cm). Pad size is 4.5"x4" (11x10cm).














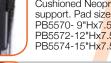


















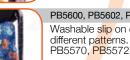



PA3062 Pelvic Positioning Strap PB3052 Pelvic Positioning Strap 2" (5cm) wide belt has a Velcro® belt with a D-Ring closure, provides hip stability for user. Fits circumference up to PA3011 45" (114cm) and PB3052 57" (145cm).



PB5592 Pelvic Padded Positioning Strap PB5594 Pelvic Padded Positioning Strap 6"x10" (15x25cm) pad attached to a 2" (5cm) Velcro® belt with a D-Ring closure. Fits circumference up to PB5592 45" (114cm) and PB5594 57" (145cm).



PB5588, PB5590 Lateral Supports Keeps the upper torso in symmetrical

alignment. Adjustable from PB5588 7-11" (18-28cm) and PB5590 10-16" (25-41cm) with a 3.5" (9cm) height range. Pad size is 4"Hx6"W (10x15cm). For planar pads only.



PNG30170 X-Style Chest Vest-11"Lx9.5"W PNG30056 X-Style Chest Vest-14.5"Lx9.5"W Neoprene® with camlock adjustment provides anterior trunk support for child.

PNG30171 Y-Style Chest Vest-17.5"Lx11.5"W

PNG30021 Y-Style Chest Vest-20.5"Lx12.5"W

Neoprene® with camlock adjustment

provides anterior trunk support for child.

PB5596 Upper Body Padded Positioning Strap

PB5598 Upper Body Padded Positioning Strap

6"x10" (15x25cm) pad attached to a 2"

(5cm) Velcro® belt with a D-Ring





PB5554 Black Molded Tray-21"Lx24"W Trav is accessible in all positions and provides anterior support as user moves from supine to upright and prone positions. Adjustable in depth and height.

PB5556 Black Molded Swing-Away Tray-21"Lx24"W

Swing-away tray is accessible in all positions and provides anterior support as user moves from supine to upright and prone positions. Adjustable in depth, height and tray angle. Tray size is 21"Lx24"W (53x61cm).

PB5558 Clear Swing Away Tray-21"Lx24"W

Swing-away tray is accessible in all positions and provides anterior support as user moves from supine to upright and prone positions. Adjustable in depth, height and tray angle. Tray size is 21"Lx24"W (53x61cm).

PT50156 Elbow Pad-9.25" Cutout

Additional padded support surface on the black molded or clear tray. 9.25" (24cm) cutout width. Overall width 21.5" (55cm), pad extends back 3.75" (10cm) posterior from current table edge for additional upper body/extremity support.

PT50160 Elbow Pad-11" Cutout

Additional padded support surface on the black molded or clear tray. 11" (28cm) cutout width. Overall width 21.5" (55cm), pad extends back 3.75" (10cm) posterior from current table edge for additional upper body/extremity support.

PY5582 Elbow Pad-13" Cutout

Additional padded support surface on the black molded or clear trav. 13" (33cm) cutout width. Overall width 21.5" (55cm), pad extends back 3.75" (10cm) posterior from current table edge for additional upper body/extremity support.

PNG30295 Hygienic Padded Cover

Washable slip on cover is padded and available in five different patterns.

PNG50042 Hand Grips

Provides additional arm and hand support and alignment. Easily movable for precise arm positioning. Attaches with suction cups. Clear tray required.

PB5560 Push Handle Makes pushing the stander easier for the caregiver and also aids in lifting and lowering the child. Required for dual control.

NEW TANK

P80809 Tool Pouch Holds the assembly tools provided with the stander. Attaches to unit with Velcro® straps.

standing progress and leg abduction.

P80246 Angle Locator Magnetic angle locater can be utilized for





PB3049 Hygienic Cover Washable slip on cover available in five different patterns. For PB5588 and PB5599.




PB5582 Head Support Planar Pad Cushioned Neoprene® pad for head support. Pad size 9"Hx10"W (23x25cm).



PY5628 Head Support 6"Hx10"W PY5630 Head Support 7"Hx14"W Cushion sizes are: PY5626-5"Hx8"W (13x20cm), PY5628-6"Hx10"W (15x25cm), PY5630-7"Hx114"W (18x36cm). Depth range is +1" to -2.75" (+2 to -7cm) from backrest. Maximum height from top of upper body pad is 10" (25cm).



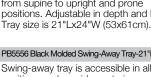

PA5586 Head Support-Form to Fit

upper body pad is 10" (25cm).

different patterns. For PB5582.

Planar Pad

and PY5630.

PA3037 Hygienic Cover

Cushion size 5"Hx10"W (13x25cm), Depth

range is -2.5" (6cm) to +1.5" (4cm) from

PB3038 Hygienic Cover for Head Support

Washable slip on cover available in five

PA3035, PA3036, PB3033 Hygienic Covers

Washable slip on cover available in five

different patterns. For PY5626, PY5628

Washable slip on cover available in five

different patterns. For PA5586.

backrest. Maximum height from top of





page 20

Zing Prone



The Prone base includes:

- Four Swivel Locking Casters
- Foot Operated Gas Spring or Pow'r Up Lift
- Mast Receiver



PB5606 Zing Prone Base Gas Spring Lift Comes with four swivel locking casters, foot operated gas spring, and a mast receiver.



PB5608 Zing Prone Base Pow'r Up Lift Comes with four swivel locking casters. Pow'r Up Lift, mast receiver.



PB5512 Mast Mast is the base for support and positioning components.



PB5514 Mast with Leg Abduction Mast is the base for support and positioning components and allows for 0-30° of independent leg abduction. Multi-adjustable knee pads and multi-adjustable foot plates required for leg abduction.

PB5640 Dual Control



Allows both foot and hand control of the gas spring lift. Push handle PB5632 required.



PA5614 Gas Spring Lift Lockout Locks gas spring in any position from lying down to fully upright and supine positions.



PB5510 Directional Locking Caster Great for transporting stander. Directionally locks one caster while leaving the other three free to swivel.















PB5622, PB5624, PB5626 Multi-Adjustable Knee Pads

Multi-adjustable knee pads are adjustable in swivel, height +/-.75" (2cm), PA5622 inside width 3.25" (8cm), PA5624 inside width 4.25" (11cm) and PB5626 inside width 5" (13cm). With a depth range of 7" (18cm).



PB5636 Knee Strap

4"x6" (10x15cm) pad attached to a 1.5" (4cm) Velcro® belt with a D-Ring closure. Fits circumference up to 19" (48cm). Adjustable in height. For PB5622, PB5624, and PB5626.

Size 2 -

154 lbs (70 kg)

PB3000 Planar Pad-10"Hx11"W

PB3002 Planar Pad-10"Hx13"W Cushioned pad covered with Neoprene®. Pad sizes: PB3000 10"Hx11" (25x28cm) PB3002 10"Hx13"W (25x33cm).



PB3004 Form to Fit Pad-10"Hx11"W PB3006 Form to Fit Pad-10"Hx13"W

Form to fit pad can be custom fit to the user. Pad sizes: PB3004 10"Hx11"W (25x28cm) PB3006 10"Hx13"W (25x33cm).


PB5538, PB5540 Hip Supports

Keeps the hips in symmetrical alignment. Adjustable from PB5538 7"-13" (18-33cm) and PB5540 12"-18" (30-46cm) with fore/aft 4" (10cm). Pad size is 4.5"x4" (11x10cm). For planar pads only.

PA3011, PT30066 Pelvic Positioning Strap

2" (5cm) wide belt has a Velcro® belt with a D-Ring closure, provides hip stability for user. Fits circumference up to PA3011 29" (74cm) and PT30066 35" (89cm).

PA5610, PB5542

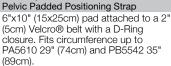

PB3025, PB3026 Hygienic Cover

Washable slip on cover available in five different patterns. For planar pads PB3000 and PB3002.





PB5612 Platform Foot Plate

positioning of the feet.

PB5616, PB5618, PB5620

PA5542, PA5544, PB5526

Multi-Adjustable Foot Plates

straps can be added for more precise

Foot plates are adjustable forward/aft +/-

and in height. Foot holder sizes: PB5616 7.75"Lx3.25"W (20x8cm), PB5618 9.75"Lx4"W (25x10cm). PB5620 11.75"Lx5"W (30x13cm).

Foot Holders for Platform Foot Plate

.75" (2cm) +extra 1.5" (3.8cm), plantar/dorsi

+/- 20°, toe-in/toe-out, swivel 30° left/right,

PY5638 Secure Foot Straps-15"L PNG50506 Secure Foot Straps-21"L Velcro® and D-ring adjustment holds feet securely in place at ankles and toes with multiple attachment slots for proper foot positioning. Straps measured in length over top of foot, with approximately 1.5" (4cm) of Velcro® overlap.

PT50080 Foot Straps-10"L PY5636 Foot Straps-15"L PNG50505 Foot Straps-21"L

Velcro® and D-ring adjustment holds feet in alignment with multiple attachment slots for proper foot positioning. Straps measured in length over top of foot, with approximately 1.5" (4cm) of Velcro®

PB5528 Calf/Knee Pad with Strap

17.5"x5.5" (44x14cm) pad attached to a 2" (5cm) Velcro® belt with a D-Ring closure, provides calf/knee stability for user. Fits circumference up to 28"-38" (71-97cm). Gives support from the back of the calf/ knee. Adjustable in height and fore/aft.





W



PB3027, PB3028 Hygienic Cover Washable slip on cover available in five different patterns. For form to fit pads PB3004 and PB3006.



PA3006, PA3007, PA3041, PA3042 Planar Pad Cushioned pad covered with Neoprene®. Pad

SIZES: PA3006-5"Hx7"W (13x18cm) PA3007-5"Hx9"W (13x21cm) PA3041-7"Hx7"W (18x18cm) PA3042-7"Hx9"W (18x21cm)



PA3079, PB3046, PA3081, PB3047 Form to Fit Pad

Form to fit pad with extra sternum positioning. Built in lateral support that can be custom fit to the user. Pad sizes are: PA3079-5"Hx12"W (13x30cm) PB3046-6"Hx13"W (15x33cm PA3081-7"Hx12"W (18x30cm) PB3047-10"Hx13"W (25x33cm

PB5627 Lateral Supports

Keeps the upper torso in symmetrical alignment. Adjustable range of 7"-12" (18-30cm) with 7"W planar pads, 9"-14" (23-36cm) with 9"W planar pads. For planar pads only.



PA5610, PB5542 Upper Body Padded Positioning Strap

6"x10" (15x25cm) pad attached to a 2" (5cm) Velcro® belt with a D-Ring closure. PA5610 fits circumference up to 29" (74cm), PB5542 fits circumference up to 35" (89cm).



PA3028, PA3029. PA3045, PA3046 Hygienic Cover

Washable slip on cover available in five different patterns. For planar pads PA3006, PA3007, PA3081, and PA3042,



PA3086, PB3050, PA3088, PB3051 Hygienic Cover Washable slip on cover available in five different patterns. For form to fit pads PA3079, PB3046, PA3081, and PB3047.



PB3049 Hygienic Cover for Lateral Supports Washable slip on cover available in five different patterns. For lateral supports PA5627.



Pad size is 9"Hx3"W (23x8cm). With a 7" (18cm) range from the upper body support, individual face aperture pads are



PY5582 Elbow Pad-13" Cutout

upper body/extremity support.

PT30093 Hygienic Padded Cover

Washable slip on cover is padded and available in five different patterns.

PNG50042 Hand Grips



Provides additional arm and hand support and alignment. Easily movable for precise arm positioning. Attaches with suction cups.

PB5632 Push Handle for Tray

Makes pushing the stander easier for the caregiver and also aids in lifting and lowering the child. Required for dual control.





P80246 Angle Locator

P80809 Tool Pouch

Magnetic angle locater can be utilized for standing progress and leg abduction.

Holds the assembly tools provided with the stander. Attaches to unit with Velcro®

PT50156 Elbow Pad-9.25" Cutout

PA5710 Hygienic Cover

PB5638.

Washable slip on cover available in five

PB5628 Black Molded Tray-21"Lx24"W

Adjustable in depth and height. Tray size

Tray provides anterior support. Adjustable

in depth, height and tray angle. Tray size is

Tray provides anterior support.

PB5630 Clear Tray-21"Lx24"W

21"Lx24"W (53x61cm).

is 21"Lx24"W (53x61cm).

different patterns. For face aperture

Additional padded support surface on the black molded or clear tray. 9.25" (24cm) cutout width. Overall width 21.5" (55cm), pad extends back 3.75" (10cm) posterior from current table edge for additional upper body/extremity support.

PT50160 Elbow Pad-11" Cutout

Additional padded support surface on the black molded or clear tray. 11" (28cm) cutout width. Overall width 21.5" (55cm). pad extends back 3.75" (10cm) posterior from current table edge for additional upper body/extremity support.

Additional padded support surface on the black molded or clear tray. 13" (33cm) cutout width. Overall width 21.5" (55cm). pad extends back 3.75" (10cm) posterior from current table edge for additional

page 22



straps.









EasyStand

Altimate Medical, Inc. 262 West First Street Morton, MN 56270 USA Tel: 800.342.8968 • 507.697.6393 Fax: 877.342.8968 • 507.697.6900 info@easystand.com • easystand.com

Standing technology should only be used under the guidance of a physician with recommendations for standing program protocol and any medical precautions. Standing programs should be monitored by the attending therapist. AMI maintains a policy of continual product improvement and reserves the right to change features, specifications, and prices without prior notification. Check with AMI for latest information. FORM ZINGCAT050918 REV C Copyright © 2018 Altimate Medical, Inc. All rights reserved. Printed in the USA.