

# CS8 single floorline bed

PERFECTLY STYLED



The Invacare CS8 bed is an electrically actuated care bed with a range of functions and positional movements that assist caregivers and provide for patient comfort. Ease of use, functionality, serviceability and cleaning practicality are key features of the bed.



# CS8 floorline bed

The bed has a fabricated steel structure with a powder coated finish. The plastic components are made from urethane, PP, POM, nylon and ABS plastic.

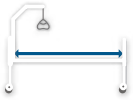
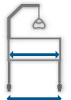
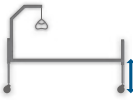


## Clinical Applications and Area of Use:

The bed is appropriate for adult patient care under supervised medical monitoring and provides support with various positioning options.

## Features

- Total safe working load (SWL) of the bed 250Kg.
- 75mm twin Tente castors (four braking castors and four free-wheeling).
- Electrically adjusted height range. The bed works within a height range of 150mm low height to a maximum high height of 680mm.
- Trendelenberg head and foot down tilt 14°.
- Four sectioned electronically operated mattress platform with independent adjustment of backrest angle (maximum 65°) and knee section angle (maximum 35°).
- Four soft contoured mattress retainer handles.
- Integral foot end mattress base extension. Extends the bed by 195mm.
- (Optional) bolt on 3/4 length folding side rails.
- Removable timber head and foot panels.
- One button chair positioning.
- Under bed lighting.

## Tech Specs

					
	Length	Width	Adjustable Height	User weight	Bed weight
CS8 floorline Single	2110mm	920mm	150 - 680mm	250kg	119kg
Mattress*	1970mm	880mm			

### FOR FURTHER INFORMATION:

**New Zealand:** Freephone 0800 468 222, email [sales@invacare.co.nz](mailto:sales@invacare.co.nz) or visit [www.invacare.co.nz](http://www.invacare.co.nz)

**Australia:** Freephone 1800 460 460, email [orders@invacare.com.au](mailto:orders@invacare.com.au) or visit [www.invacare.com.au](http://www.invacare.com.au)

All rights reserved. All information quoted is believed to be correct at time of publication. Invacare reserves the right to alter product specifications without prior consent. 04/2024



Yes, you can.®