

ROYALE

Pro Shift Patient Mover

Sit-to-Stand Transfer Aid

Effective support for promoting mobilisation, independence and patient transfer comfort.



A supportive standing aid that promotes user independence while enabling single-caregiver transfers. Ideal for weight-bearing users who need extra assistance when transferring between seating surfaces. Features ergonomic design elements for both user and carer safety.


Pro Shift Patient Mover

Why Choose the Royale Pro Shift Patient Mover?

- **Enhanced Maneuverability:** Lightweight base design with smooth 10cm castors enables effortless positioning and transport in any care environment.
- **Ergonomic Control:** Dedicated attendee push bar and adjustable foot pedals for leg spread provide caregivers precise control around furniture and fixtures.
- **User Comfort:** Padded support bar and contoured swing-away seat pads ensure comfortable, dignified transfers.
- **Adaptable Design:** Height and angle adjustable knee pad with flexible leg spread accommodates various user needs and environments.
- **Versatile Access:** Adjustable base width enables easy positioning around toilets, commodes, and chairs for seamless transfers.



Padded Lower Knee Support

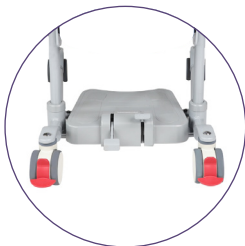

Locking Rear Castors
Ergonomic Base Design with Low Profile



Comfortable Swing-Away Seat Pads

Specifications

	S332019
Overall Width (External)	66cm - 94cm
Overall Width (Internal)	50cm - 78cm
Overall Depth	100cm
Overall Height	104cm
Base Plate Height	10cm
Rear Width	55cm
Knee Pad Measurement (W x H)	43cm x 25cm
Knee Pad Height from Base Plate	20cm - 24cm
Handle Height from Base Plate	90cm
Seat Height from Base Plate (B x T)	56cm x 68cm
Seat Width	58cm
Seat Depth	31cm
Front Castors	10cm
Rear Castors (w/brakes)	10cm
Underbed Clearance	15cm
Attendee Push Handle Height from Floor	104cm
Safe Working Load (SWL)	200kg