

2701 Falls Monitor with SafeTime Support

cura¹
Helping To Care



***... when infection control is critical,
SafeTime prevents expiry over-runs***

2701 Falls Monitor with SafeTime Support

SafeTime bed & chair pads have an internal electronic timer which causes a "Change Pad" notification on compatible monitors, when the pad reaches its pre-programmed functional life.

SafeTime Bed & Chair Pads

- Available with 90 day expiry dates.
- Supplied with a 2480 USB breakcord.
- Made from non-allergenic soft PVC.
- Vent holes in the neck reduce the risk of moisture ingress.
- Easy to clean.



2622 SafeTime Bed Pad - 90 Days

With pressure sensors and internal timer built in, it is placed under the bed sheet to alert caregivers when a patient exits the bed.



2620 SafeTime Wheelchair / Chair Pad - 90 Days

alert caregivers when a patient exits a chair.

WHAT'S IN THE BOX (2701)



2701 Falls Monitor



2742 Soft Monitor Cover



2642 Caregiver Reset Key




Batteries



QuickStart User Guide

2701 Falls Monitor with SafeTime

- Tamper-proof reset option.
- Connects with nurse call system.
- Operates on batteries or optional external plug pack (2406).
- Works with Premium (corded) and SafeTime pads & mats.
- Low battery & input lost notification (pad lost, nurse call, AC lost).
- Change pad (SafeTime) indication.

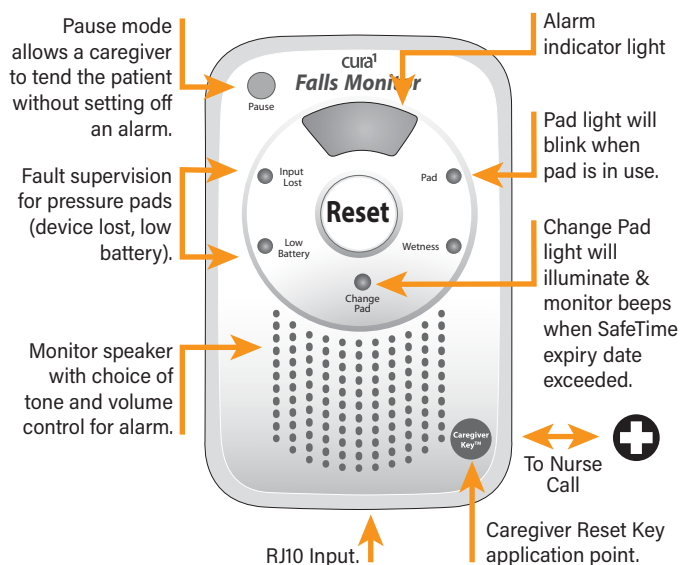

COMPATIBLE PRODUCTS

- All cura¹ premium pads, mats with RJ10 plugs.
- All cura¹ universal pads & mats.
- 2622, 2630 Safetime Pads.
- 2406 Plug Pack 6 V.



Certifications may vary between products.

cura¹
Helping To Care

10 Gray Street, Kilkenny, South Australia 5009

sales@cura1.com | www.cura1.com

+61 (0) 8 8268 8343 (International)

1300 12-CARE (2273)