

QUANTUM® R-TRAK

REAR-WHEEL DRIVE POWER CHAIR

R-TRAK

shown in *Lemon Crush* with
TRU-Balance® 3 power positioning, Q-Logic
3 Advanced Drive Control System and
Stealth® Premiere cushion.



AVAILABLE WITH 5.6 KPH
10"
SAFE SEAT ELEVATION



NEXT LEVEL REAR-WHEEL DRIVE

The R-TRAK redefines rear-wheel drive power chair performance, delivering a smooth, responsive ride with independent suspension on all four wheels. Designed for those who prefer the power, control, and natural tracking of rear-wheel drive, its compact base offers improved manoeuvrability without compromising stability.

Standard features include fender lights, a USB charger, and drive wheel colour accents. Optional 10" seat elevation allows users to drive up to 5.6 km/h while elevated—enhancing access, visibility, and confidence in motion.


QUANTUM®



| quantumrehab.com.au

Specifications

DRIVE WHEELS	35.6 (14") Pneumatic or Solid Tyres 14x4" wide Flat Free Tyres
CASTER WHEELS	Front: 22.8 cm (9") solid or Pneumatic
ANTI-TIP WHEELS	Rear: 10 cm (4") Solid
MAXIMUM SPEED ¹	Up to 10 km/h (5.6 km/h when elevated)
GROUND CLEARANCE	10 cm (3.97") (battery tray)
TURNING RADIUS ²	65.2 cm (25.7")
OVERALL LENGTH ²	88.3 cm (34.8")
BASE WIDTH	61.4 cm (24.2")
SEATING SIZES	
TRU-Balance® 3 Static Seat	Width: 12-24" (30.5 cm - 61 cm) Depth: 12-22" (30.5 cm - 55.8 cm)
TRU-Balance® 3 Power Tilt	Width: 12-24" (30.5 cm - 61 cm) Depth: 12-22" (30.5 cm - 55.8 cm)
TRU-Balance® 3 Power Lift & Tilt	Width: 12-24" (30.5 cm - 61 cm) Depth: 12-22" (30.5 cm - 55.8 cm)
SEAT-TO-FLOOR HEIGHTS	
TRU-Balance® 4	N/A
TRU-Balance® 3 Static Seat	43.6 cm - 48.7 cm (16 - 18")
TRU-Balance® 3 Power Tilt ⁵	43.6 cm - 48.7 cm (17.2 - 19.2")
TRU-Balance® 3 Tilt & iLevel® ⁵	43.6 cm - 48.7 cm (17.2 - 19.2")
BATTERY SIZE ³	MK 70Ah(2)
BATTERY WEIGHT ⁷	23.3 kg each
AVAILABLE ELECTRONICS ⁸	120A Q-Logic 3
RANGE PER CHARGE ³	Up to 26.3 km
BATTERY CHARGER	8A Off-board
MOTOR PACKAGES	4-pole motors
WEIGHT CAPACITY ⁴	150 kg
BASE WEIGHT	201 kg

- 1) Speed varies with user weight, terrain type, battery charge, battery condition and tire pressure.
- 2) Without front riggings
- 3) Battery capacity may change when adding power positioning systems.
- 4) Weight capacity is dictated by the selection of seating.
- 5) All measurements +/- 3 cm with 45.7 cm x 45.7 cm seat, fixed towers, 136 kg. weight capacity and solid tires (add 3/8" with pneumatic tires). Widths may increase seat-to-floor height. Contact Customer Service for details. Listed measurements are in the lowest position; does not include elevating system unless specified.
- 6) Battery weight may vary +/- 1.3 kg based on manufacturer.
- 7) All Quantum Rehab® units can be ordered with specialty controls.
- 8) Call for availability
- 9) Seat-to-floor heights are measured with the unit unloaded. Due to the R-Trak's suspension system, seat-to-floor heights are approximately 0.63 cm lower when the unit is loaded.


Standard Features

- 4-pole motors
- Maximum speed of up to 10 km/h (5.6 km/h when elevated)
- LED lighting in front caster arms and fenders
- Superior, automotive-grade SRS (Smooth Ride Suspension)
- Unmatched driving performance and ride quality
- Advanced stability for smooth obstacle transitions
- ISO 7176-19 transit crash tested approved with all configurations

Options

- TRU-Balance® 3 Power Positioning: 10" of power seat elevation while at walking speed (5.6 km/h)
- Power Tilt, Recline, and Elevating Legrests
- USB Charger
- LED Lighting Package
- Quantum® backup camera
- 10cm (14") wide flat-free or pneumatic knobby extended traction tyres^{9,4}

Colours

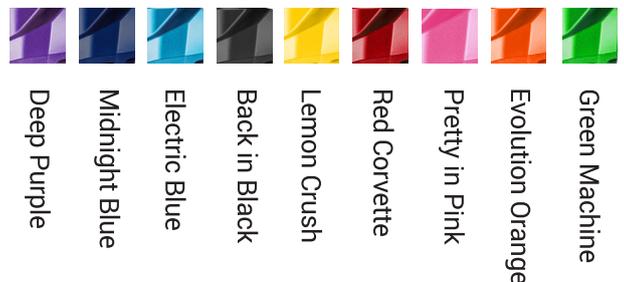

QUANTUM®
#1 FOR REHAB POWER

 **PROUDLY BUILT IN AUSTRALIA**

(03) 8770 9600 | quantumrehab.com.au